

HCP_Coverage_example_using_SpotMap_v4.0_software

Summary

Software Version: 4.0.027

Author: Michael Hutchins

Date: 10/05/2018

Example analysis report based on images donated by University of Cambridge.


Coverage of 14% determined using SpotMap v4.0 software (TotalLab Ltd).

Analysis report includes Audit Trail to support 21CFR Part 11 compliance.

Coverage: 14%

Total Number of Spots: 862


Number of Common Spots: 116

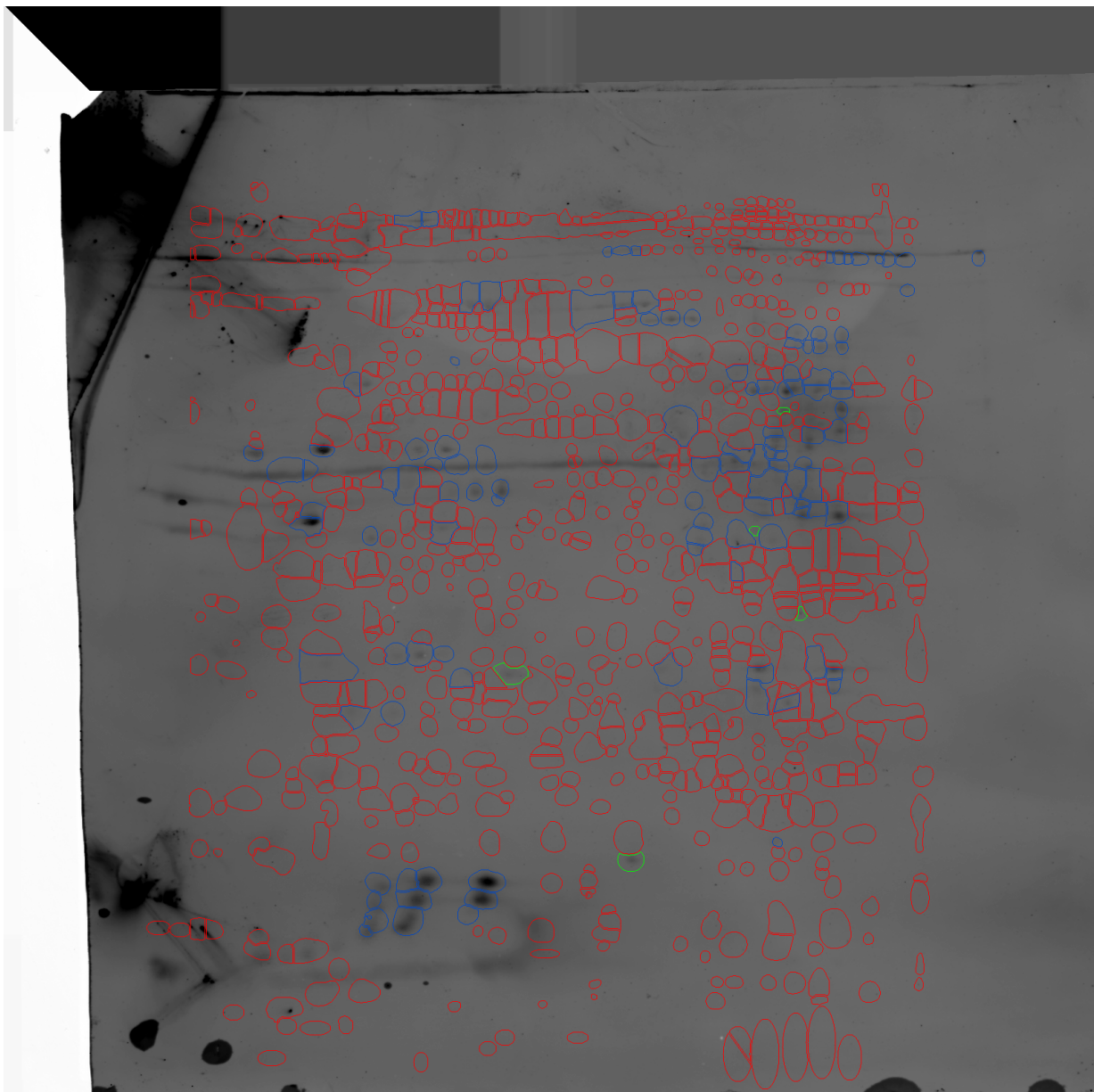

Gel Image

Gel image donated by University of Cambridge


- 1 Gel
- 1600x1600 8bit
- Nickname: n/a
- QC - - Color depth is dangerously low - 8 bits per pixel (256 unique values) makes quantification very inaccurate.
- QC - - Image is saturated.
- Present Spots: 857
- Absent Spots: 5

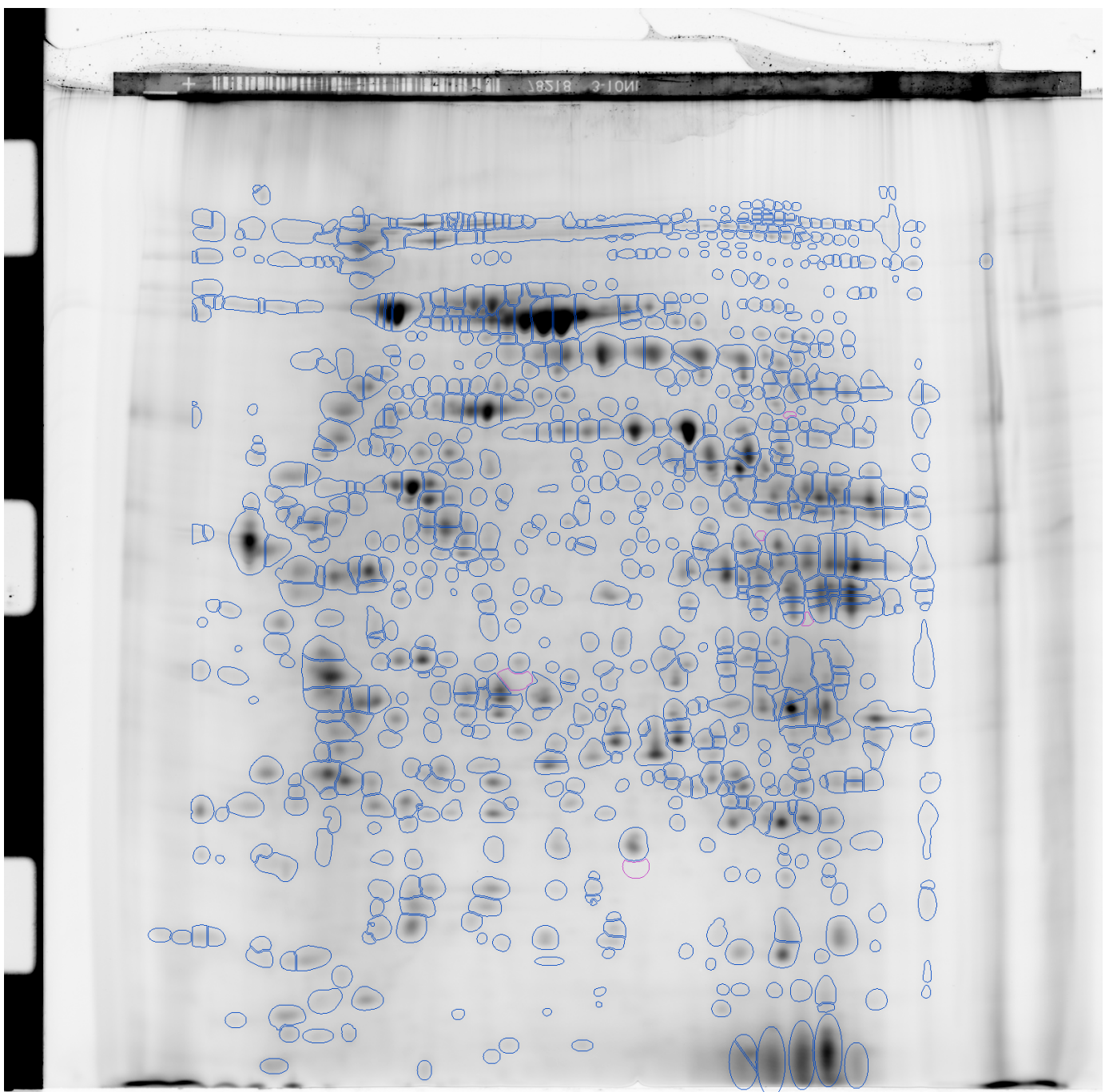

Blot Image


Blot image donated by University of Cambridge

- 2 Blot
- 1600x1600 8bit
- Nickname: n/a
- QC - - Color depth is dangerously low - 8 bits per pixel (256 unique values) makes quantification very inaccurate.
- Present Spots: 121
- Absent Spots: 741
- Number of Alignment Vectors: 17

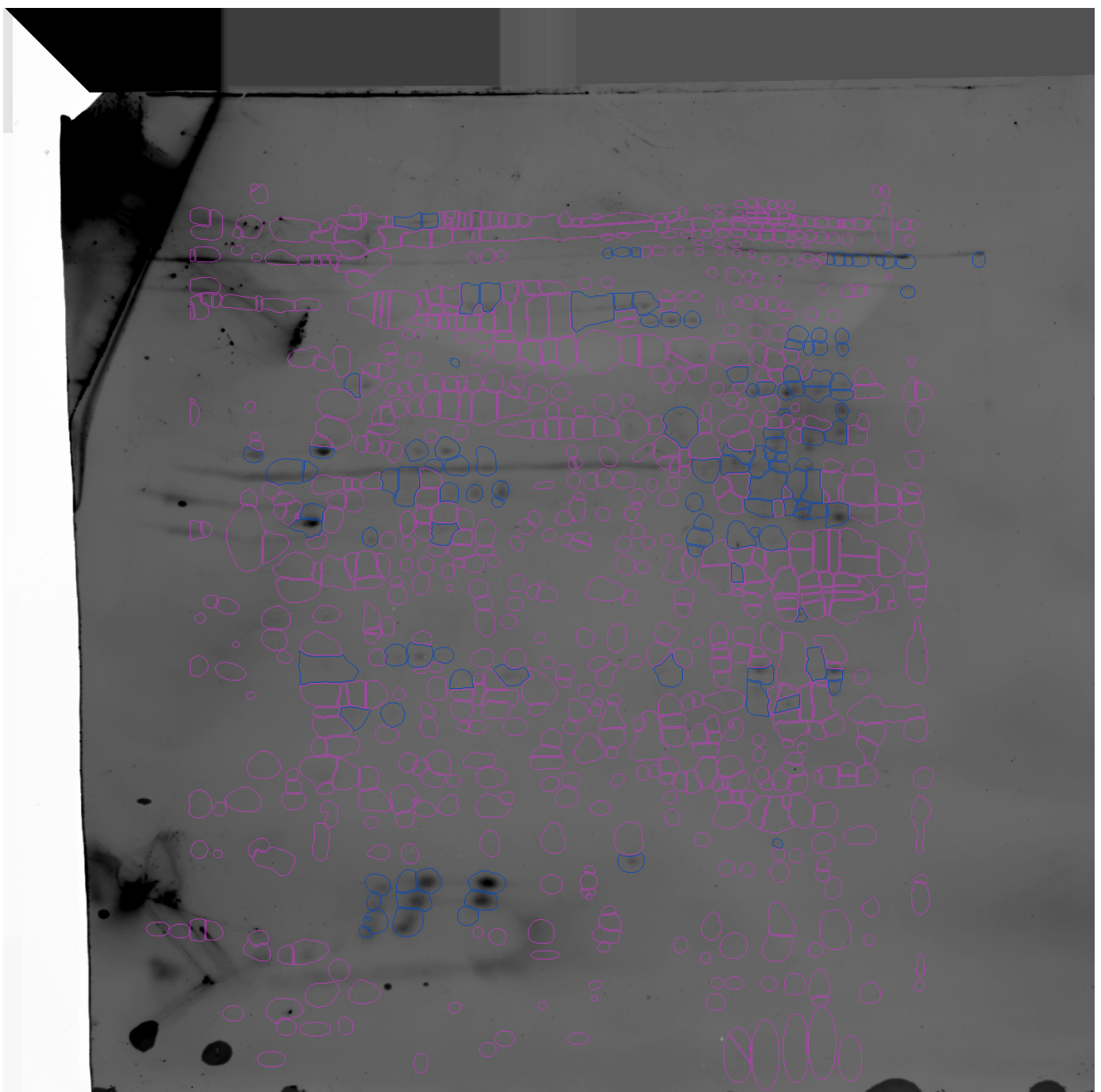


Spot Detection

- Detection Image: A: 1 Gel
- Region(s) of Interest:
 - ROI Rectangle: 272,258 1110x1232
- Peak Sensitivity: 263
- Smoothing: 1
- Original Number of Spots Detected: 795
- Final Number of Spots (inc. Edits): 862

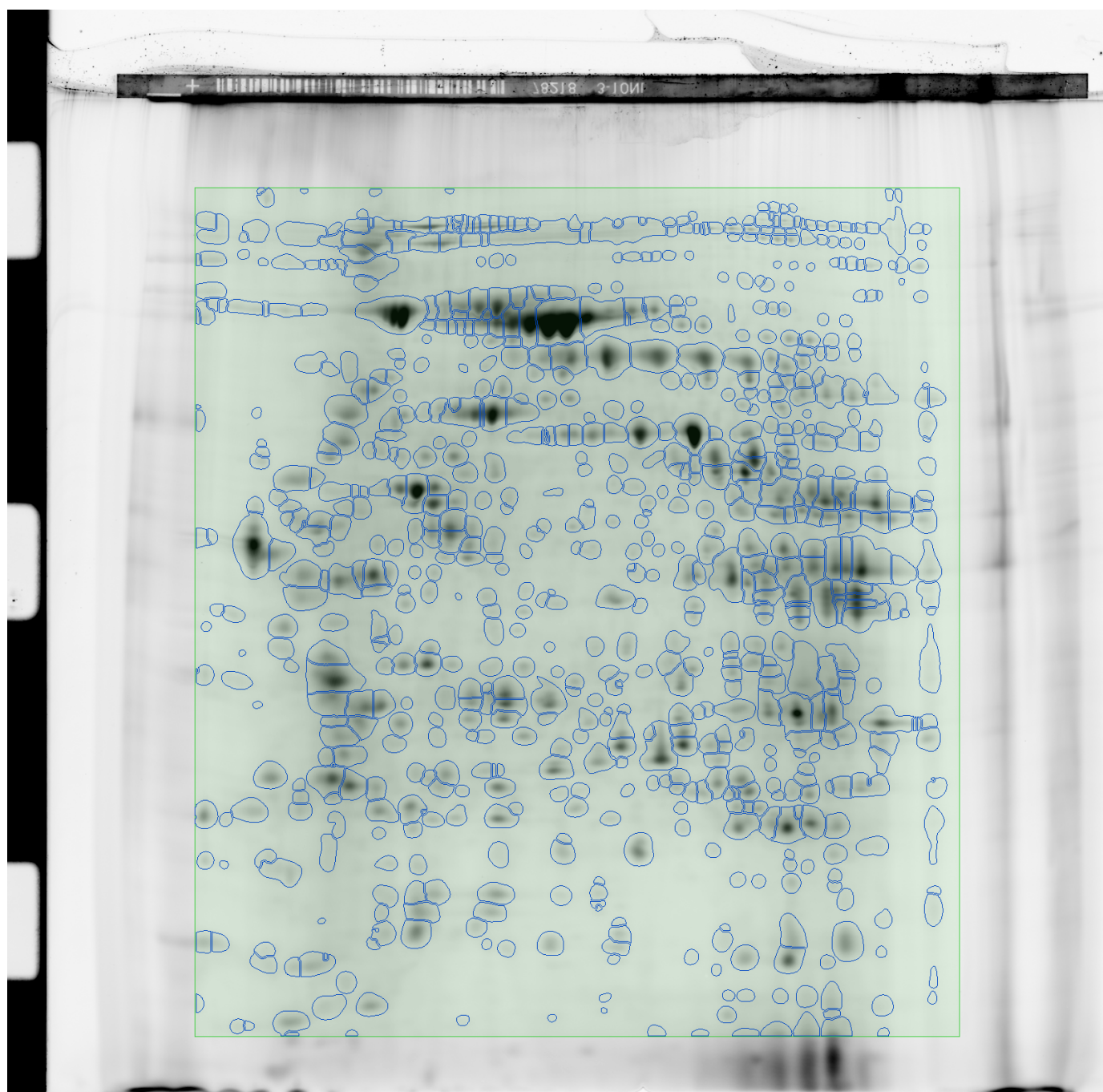


Coverage - 14%

- Number of Common Spots: 116
- Number of Unique Spots in Gel: 741
- Number of Unique Spots in Blot: 5

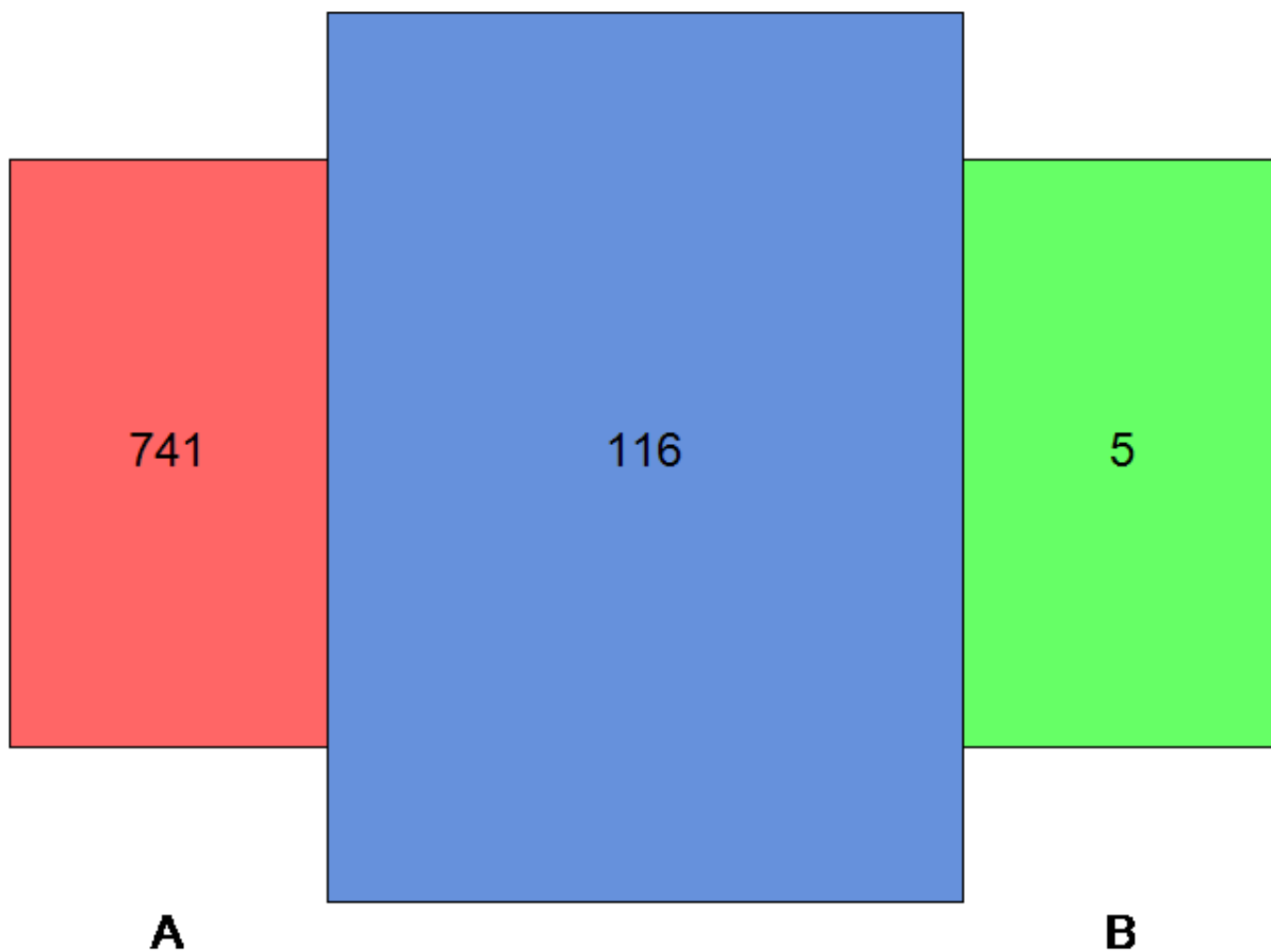


Combined Image

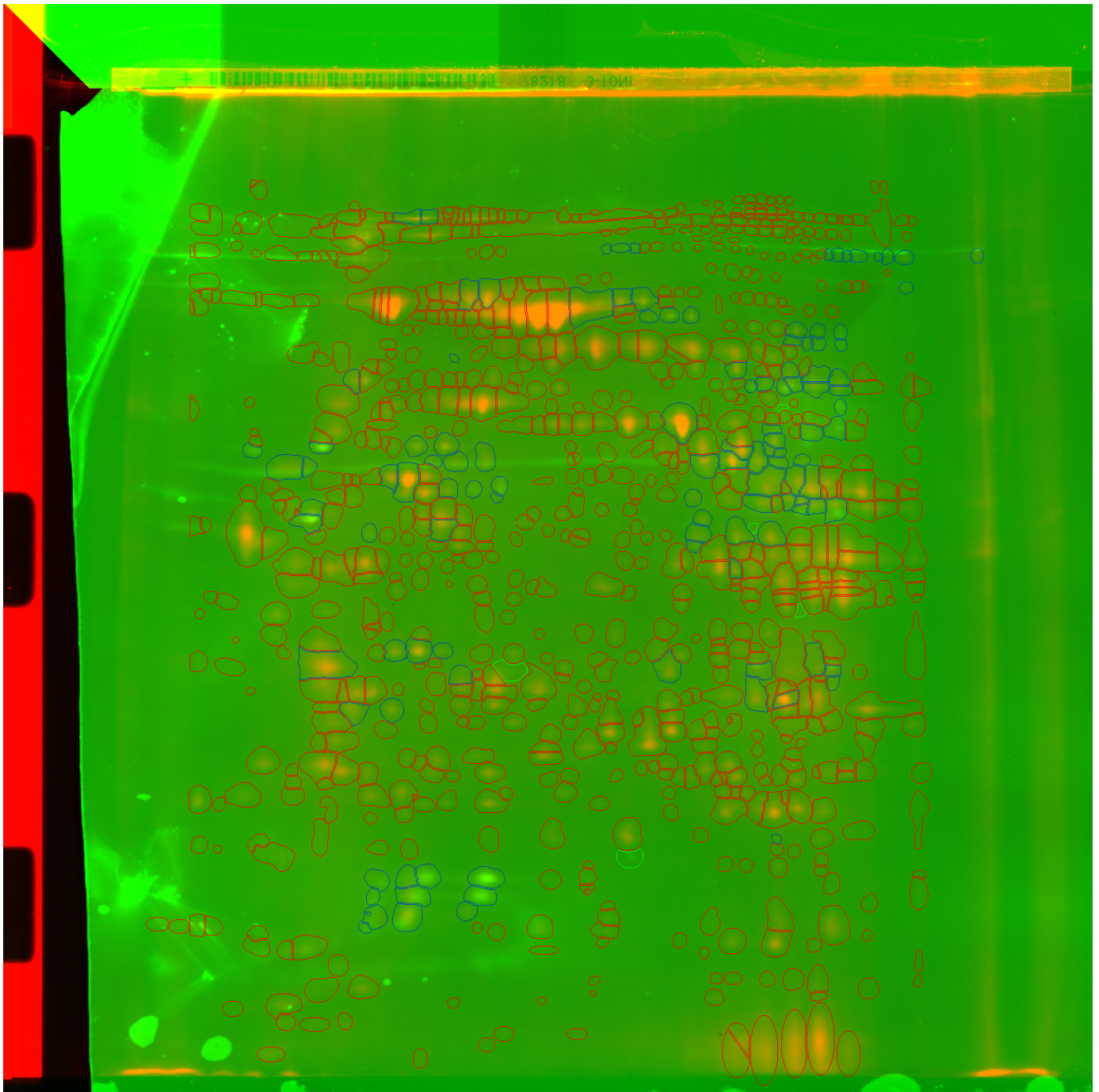


Appendix - Q1.png

B: 2 Blot

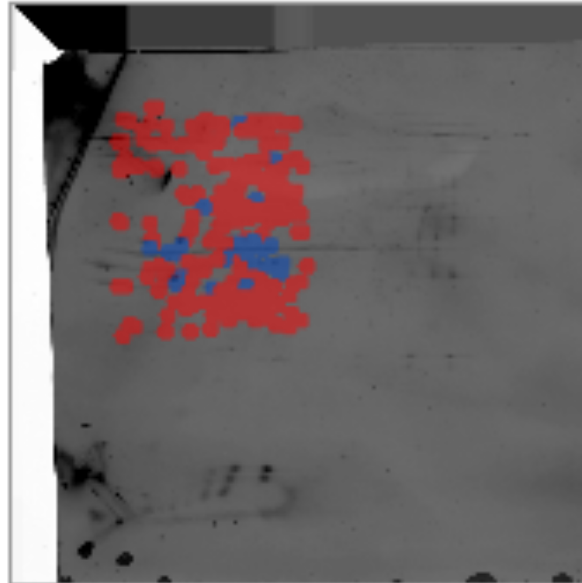


11%
coverage
■ 25 ■ 0

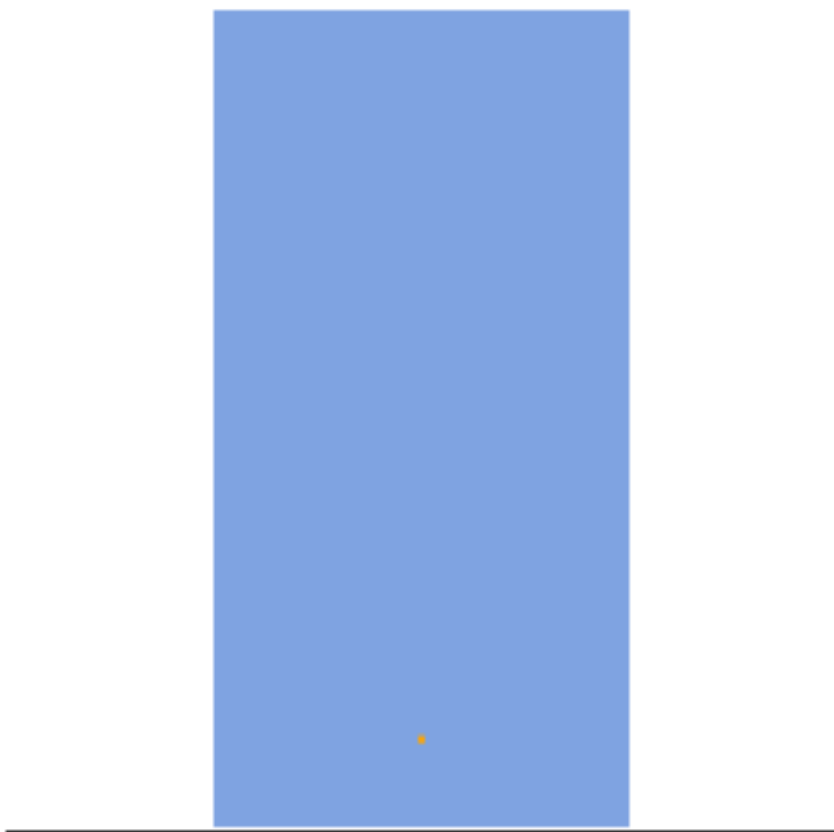


Appendix - Q2.png

B: 2 Blot

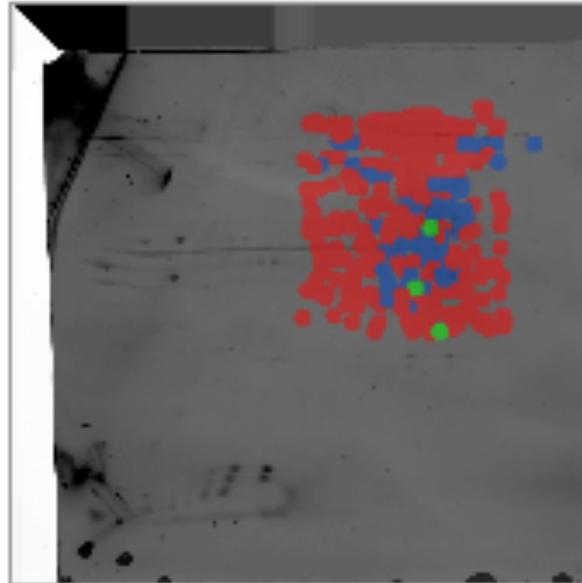


20%
coverage
■ 65 ■ 3

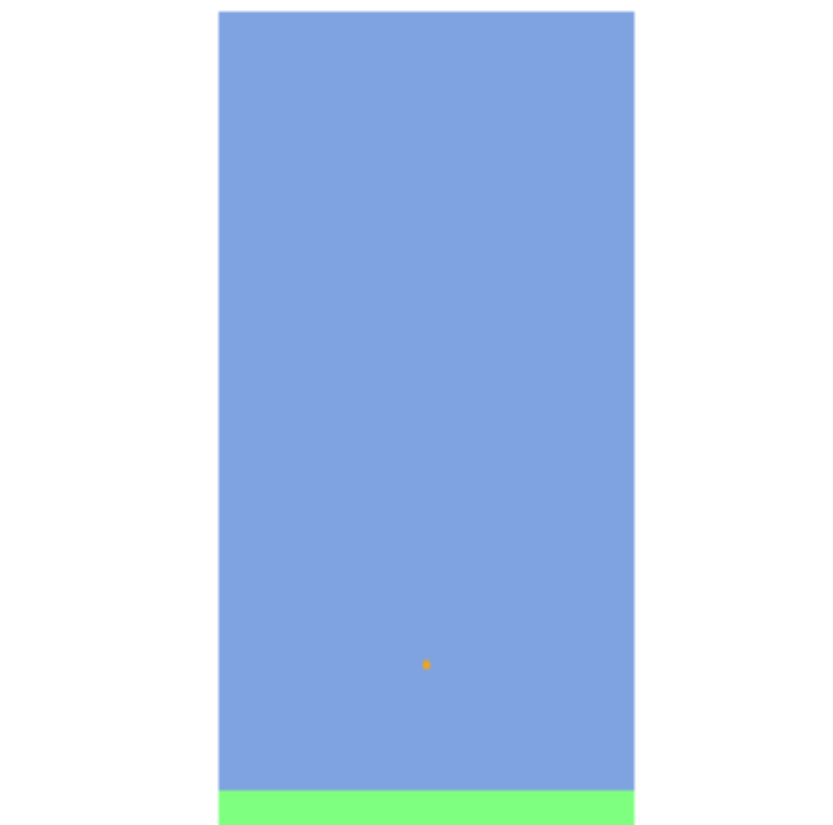


Appendix - Q3.png

B: 2 Blot

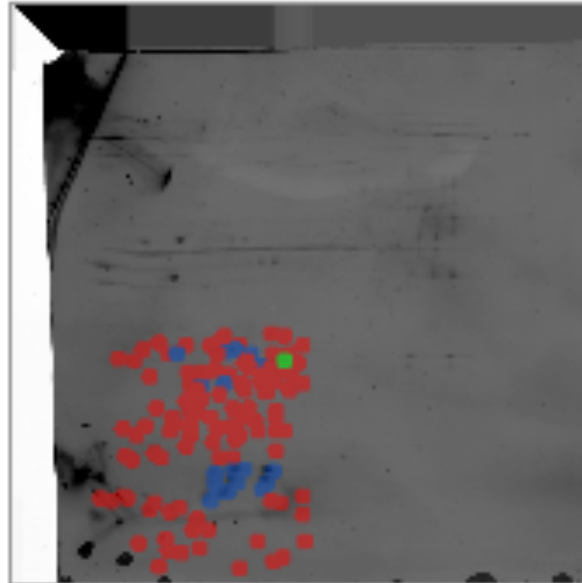


15%
coverage
■ 18 ■ 1

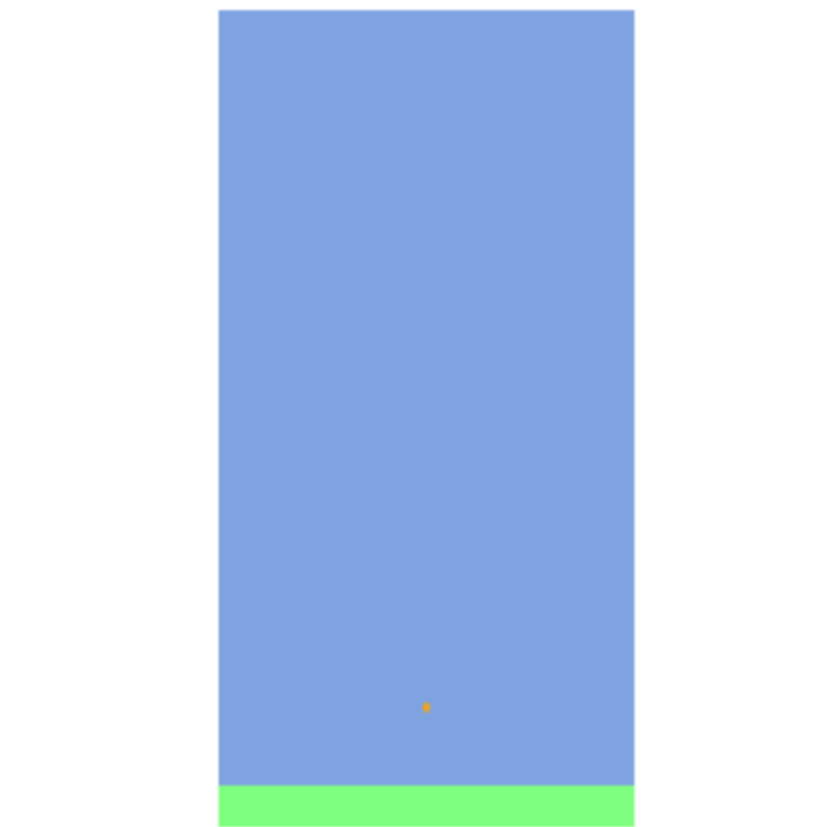


Appendix - Q4.png

B: 2 Blot

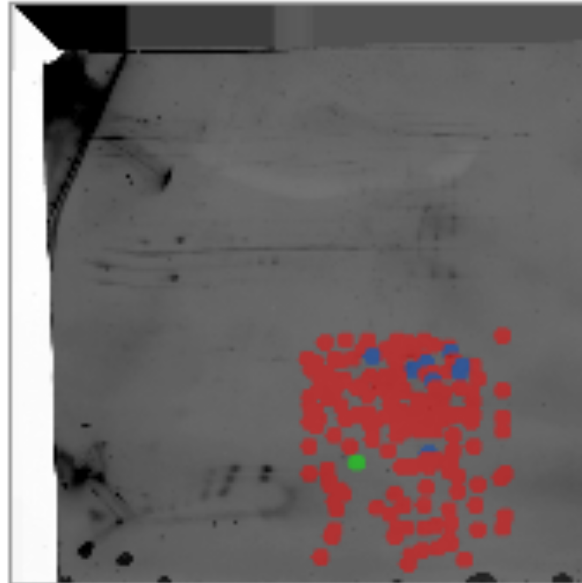

5%
coverage
■ 8 ■ 1

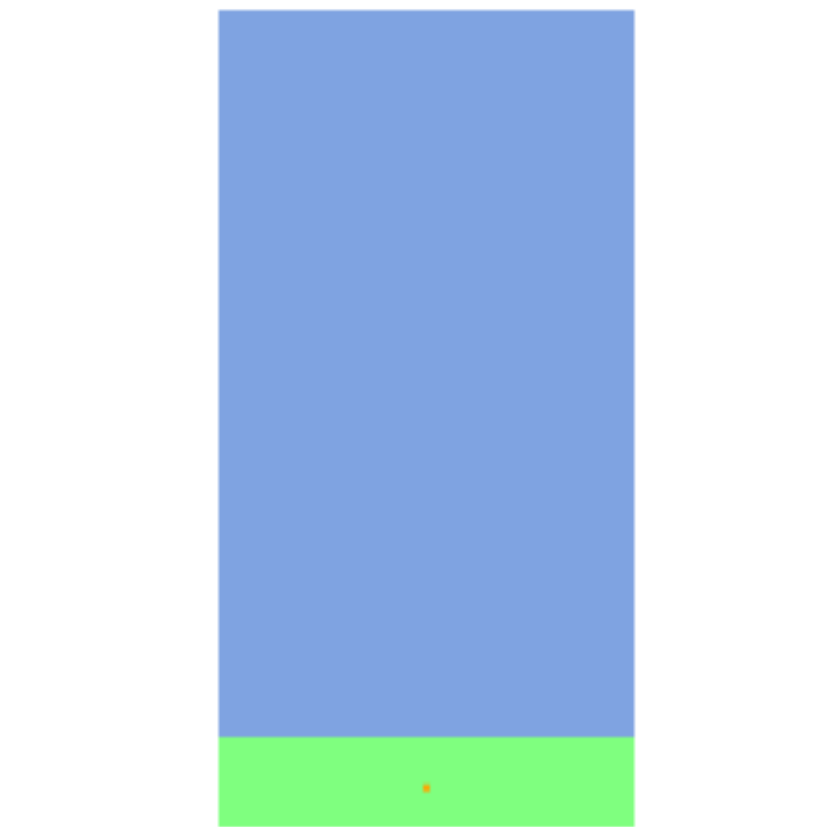

Appendix - Spot Measurements

#	Area	Present (A)	Present (B)	Volume (A)	Volume (B)	Raw Volume (A)	Raw Volume (B)	Peak Height (A)	Peak Height (B)	Peak X (A)	Peak X (B)	Peak Y (A)	Peak Y (B)	Raw Peak Height (A)	Raw Peak Height (B)	Bg Total (A)	Bg Total (B)	MW	pI
1	81	TRUE	FALSE	700	182	3697	13304	14	4	371	362	261	263	51	166	2997	13122	-	-
2	514	TRUE	FALSE	5444	1025	23948	83779	26	8	377	381	273	270	62	169	18504	82754	-	-
6	193	TRUE	FALSE	581	242	4827	29385	7	5	1280	1283	260	260	29	156	4246	29143	-	-
7	177	TRUE	FALSE	911	185	4628	26912	9	8	1292	1289	260	262	30	159	3717	26727	-	-
8	227	TRUE	FALSE	1453	201	6901	34692	15	7	1111	1106	287	292	39	159	5448	34491	-	-
9	176	TRUE	FALSE	862	159	5262	26910	12	2	1095	1096	287	279	37	154	4400	26751	-	-
10	132	TRUE	FALSE	667	131	4099	20063	13	3	1140	1141	289	285	39	154	3432	19932	-	-
11	134	TRUE	FALSE	797	171	4281	20405	14	2	1125	1125	288	282	40	153	3484	20234	-	-
12	1549	TRUE	FALSE	11265	4093	45343	234894	15	9	1282	1277	318	308	37	158	34078	230801	-	-
13	198	TRUE	FALSE	1965	384	11469	31272	18	4	519	516	300	292	66	160	9504	30888	-	-
14	679	TRUE	FALSE	5021	8026	29465	130246	13	30	295	296	310	311	49	210	24444	122220	-	-
15	162	TRUE	FALSE	821	287	6653	24911	9	6	994	999	298	293	45	158	5832	24624	-	-
19	1252	TRUE	FALSE	10731	12185	54551	228781	17	22	302	309	327	310	52	195	43820	216596	-	-
20	114	TRUE	FALSE	1286	177	5618	17619	30	3	654	653	307	304	68	156	4332	17442	-	-
21	189	TRUE	FALSE	3808	207	10990	29313	43	3	666	662	316	317	81	157	7182	29106	-	-
22	367	TRUE	FALSE	9506	1757	23819	58275	88	23	692	693	315	301	127	177	14313	56518	-	-
23	307	TRUE	FALSE	3809	402	16703	48293	25	3	824	824	316	299	67	159	12894	47891	-	-
24	282	TRUE	FALSE	7343	693	18905	43839	91	4	704	702	314	298	132	157	11562	43146	-	-
27	182	TRUE	FALSE	2934	186	10578	28032	46	3	639	641	313	312	88	156	7644	27846	-	-
28	191	TRUE	FALSE	4074	483	11714	29897	53	14	670	671	316	322	93	168	7640	29414	-	-
29	243	TRUE	FALSE	5193	272	15642	37451	71	3	715	715	314	299	114	156	10449	37179	-	-
30	213	TRUE	FALSE	1253	364	9134	32953	19	3	980	973	307	303	56	156	7881	32589	-	-
32	2218	TRUE	FALSE	70242	19796	190014	361368	91	41	532	534	341	343	145	195	119772	341572	-	-
33	706	TRUE	TRUE	20782	3288	51846	110600	118	22	608	610	315	315	162	174	31064	107312	-	-
34	188	TRUE	FALSE	4832	245	12352	29385	70	5	681	685	316	322	110	160	7520	29140	-	-
35	181	TRUE	FALSE	2751	315	10534	27827	40	3	727	727	313	300	83	155	7783	27512	-	-
37	83	TRUE	FALSE	1116	271	5764	13302	35	9	528	525	315	314	91	166	4648	13031	-	-
38	504	TRUE	TRUE	13589	1401	35261	78508	118	19	617	616	314	314	161	172	21672	77107	-	-
39	255	TRUE	FALSE	3109	653	14074	39413	27	13	755	766	312	307	70	165	10965	38760	-	-
41	196	TRUE	FALSE	1633	589	9473	30773	20	7	959	963	308	312	60	161	7840	30184	-	-
42	127	TRUE	FALSE	680	184	5125	19488	15	4	1050	1054	310	310	50	156	4445	19304	-	-
43	312	TRUE	FALSE	4199	601	17927	48025	32	11	737	742	313	315	76	163	13728	47424	-	-
45	626	TRUE	FALSE	13040	1960	41836	100867	55	10	523	509	321	321	101	168	28796	98907	-	-
48	980	TRUE	FALSE	6149	16370	43389	177090	10	91	362	354	318	315	48	255	37240	160720	-	-
49	97	TRUE	FALSE	1725	136	6187	14977	33	3	648	653	312	311	79	156	4462	14841	-	-
50	221	TRUE	FALSE	1897	388	8748	33759	27	8	1163	1168	312	305	58	159	6851	33371	-	-
51	129	TRUE	FALSE	1698	337	8664	20203	32	6	569	562	315	313	86	160	6966	19866	-	-
52	194	TRUE	FALSE	1606	403	7232	29503	30	9	1180	1175	313	306	59	159	5626	29100	-	-
53	259	TRUE	FALSE	2224	410	9217	39260	30	8	1197	1193	313	307	57	158	6993	38850	-	-
54	2119	TRUE	FALSE	16790	12775	105788	351815	16	95	462	406	322	307	58	255	88998	339040	-	-
55	271	TRUE	FALSE	1886	1039	15165	43857	15	7	479	469	322	324	64	165	13279	42818	-	-
59	213	TRUE	FALSE	1666	296	6991	32246	25	8	1216	1215	315	308	50	158	5325	31950	-	-
60	309	TRUE	FALSE	1554	573	9279	46614	11	8	1251	1255	317	308	36	157	7725	46041	-	-
62	342	TRUE	FALSE	4460	782	16088	52766	55	6	1050	1065	317	324	89	158	11628	51984	-	-
63	199	TRUE	FALSE	1307	185	6083	29836	25	8	1238	1232	317	309	49	157	4776	29651	-	-
64	278	TRUE	FALSE	1250	462	7088	41884	13	3	1325	1319	318	309	34	152	5838	41422	-	-
67	67	TRUE	FALSE	353	90	1961	10073	12	2	1227	1226	317	310	36	151	1608	9983	-	-
68	255	TRUE	FALSE	3657	557	11817	39571	56	7	1078	1089	318	324	88	160	8160	39014	-	-
69	154	TRUE	FALSE	572	227	3344	23173	12	3	1330	1337	318	312	30	152	2772	22946	-	-
70	238	TRUE	FALSE	1542	315	11300	37443	15	3	958	971	318	326	56	159	9758	37128	-	-

71	175	TRUE	FALSE	1501	639	7451	27239	22	8	1099	1095	320	324	56	160	5950	26600	-	-
72	33	TRUE	FALSE	248	25	1931	5107	19	2	656	654	317	321	70	156	1683	5082	-	-
73	153	TRUE	FALSE	1264	540	6619	23796	22	8	1118	1112	322	324	57	160	5355	23256	-	-
74	124	TRUE	FALSE	883	313	5223	19161	18	5	1137	1134	323	324	53	157	4340	18848	-	-
76	552	TRUE	FALSE	6067	1040	33115	89360	19	5	484	480	336	344	68	165	27048	88320	-	-
77	920	TRUE	FALSE	20821	2403	64061	143162	58	11	565	577	328	349	105	164	43240	140759	-	-
78	168	TRUE	FALSE	445	1106	7333	29330	5	16	341	347	330	325	46	184	6888	28224	-	-
79	763	TRUE	FALSE	20633	2503	55731	117716	90	10	632	632	339	341	136	161	35098	115213	-	-
80	225	TRUE	FALSE	2112	340	8862	35213	26	5	1075	1070	330	325	56	160	6750	34873	-	-
81	1011	TRUE	FALSE	23395	3711	70912	156372	83	10	620	587	339	347	130	161	47517	152661	-	-
82	546	TRUE	FALSE	6275	1722	31391	85260	30	20	661	685	336	326	76	173	25116	83538	-	-
83	186	TRUE	FALSE	1548	382	10661	28840	18	18	696	687	334	324	67	171	9113	28458	-	-
84	264	TRUE	FALSE	3265	760	11449	41152	41	9	1097	1092	331	326	72	162	8184	40392	-	-
85	187	TRUE	FALSE	1312	185	7296	28983	15	2	1053	1052	328	325	47	156	5984	28798	-	-
86	285	TRUE	FALSE	4641	813	13761	44133	53	9	1117	1112	333	327	85	161	9120	43320	-	-
87	324	TRUE	FALSE	4555	766	15247	50014	51	7	1136	1130	334	327	84	159	10692	49248	-	-
88	213	TRUE	FALSE	1494	712	7671	32875	15	9	1183	1175	333	328	44	160	6177	32163	-	-
89	238	TRUE	FALSE	1193	659	8333	36359	11	9	1197	1193	335	329	41	159	7140	35700	-	-
90	238	TRUE	FALSE	1774	637	7962	36337	16	10	1215	1214	337	331	42	160	6188	35700	-	-
91	228	TRUE	FALSE	1588	436	7288	34636	17	9	1237	1234	338	332	42	159	5700	34200	-	-
92	174	TRUE	FALSE	965	162	4445	26262	9	3	1326	1330	339	348	29	153	3480	26100	-	-
93	877	TRUE	FALSE	31048	3313	80160	142756	88	10	517	525	352	356	144	169	49112	139443	-	-
94	223	TRUE	FALSE	1043	1185	9294	39541	9	43	344	339	355	349	46	215	8251	38356	-	-
95	250	TRUE	FALSE	984	587	11234	38837	8	6	954	970	356	353	49	159	10250	38250	-	-
96	770	TRUE	FALSE	10584	22312	36764	160142	26	76	297	282	363	367	60	255	26180	137830	-	-
97	199	TRUE	FALSE	951	596	9110	31043	9	7	938	944	356	354	50	160	8159	30447	-	-
98	158	TRUE	FALSE	704	519	6550	24693	10	7	990	985	356	352	47	160	5846	24174	-	-
99	324	TRUE	TRUE	2138	1095	17042	50667	11	10	905	918	358	357	57	163	14904	49572	-	-
100	197	TRUE	TRUE	992	903	9463	31044	9	11	928	926	357	356	52	164	8471	30141	-	-
101	187	TRUE	TRUE	887	834	9863	29258	9	11	885	885	359	360	57	163	8976	28424	-	-
102	396	TRUE	FALSE	2115	1142	21915	60938	13	10	708	710	365	369	63	161	19800	59796	-	-
103	112	TRUE	FALSE	1203	3314	5235	24482	18	66	279	278	362	367	54	255	4032	21168	-	-
104	841	TRUE	FALSE	7300	14723	37576	152647	16	57	408	397	361	371	52	221	30276	137924	-	-
105	1859	TRUE	FALSE	39984	18327	138511	295318	44	67	528	532	374	365	97	216	98527	276991	-	-
106	208	TRUE	FALSE	1005	508	11613	31916	10	7	730	728	365	369	61	158	10608	31408	-	-
107	298	TRUE	TRUE	1910	3871	9062	48273	14	49	1228	1231	369	362	38	198	7152	44402	-	-
108	168	TRUE	FALSE	740	347	8972	25883	8	11	687	690	366	372	57	163	8232	25536	-	-
109	355	TRUE	TRUE	2595	3647	10050	56187	13	47	1296	1292	372	363	34	195	7455	52540	-	-
110	169	TRUE	TRUE	822	1452	5385	26802	10	42	1208	1216	367	361	37	192	4563	25350	-	-
111	153	TRUE	FALSE	844	1530	6505	27540	8	30	372	376	364	371	45	200	5661	26010	-	-
112	163	TRUE	FALSE	827	842	5391	25454	9	33	1206	1206	367	360	37	184	4564	24612	-	-
113	405	TRUE	TRUE	4224	4365	13944	64710	26	47	1252	1250	370	362	50	196	9720	60345	-	-
114	173	TRUE	FALSE	1518	2221	9476	30247	16	37	474	467	369	373	62	199	7958	28026	-	-
115	504	TRUE	TRUE	6192	4924	15768	80524	37	69	1323	1316	372	365	56	219	9576	75600	-	-
116	177	TRUE	FALSE	1369	2870	8626	31367	15	42	465	456	370	372	56	203	7257	28497	-	-
117	201	TRUE	FALSE	1519	2177	11971	34136	14	35	488	480	368	373	66	194	10452	31959	-	-
118	172	TRUE	TRUE	1578	1550	5706	27178	18	49	1245	1238	371	362	42	198	4128	25628	-	-
119	277	TRUE	FALSE	2080	4395	12606	49269	14	43	452	435	369	373	52	205	10526	44874	-	-
120	109	TRUE	FALSE	339	190	4045	16649	6	2	1118	1117	368	365	40	153	3706	16459	-	-
121	139	TRUE	TRUE	646	822	4121	21533	8	46	1283	1282	370	364	33	195	3475	20711	-	-
122	140	TRUE	FALSE	609	253	5229	21393	9	4	1141	1138	370	366	42	155	4620	21140	-	-
123	132	TRUE	FALSE	491	391	4715	20191	8	7	1165	1172	372	373	40	157	4224	19800	-	-
124	203	TRUE	FALSE	896	925	8001	31172	8	6	1037	1037	387	380	43	155	7105	30247	-	-
125	471	TRUE	FALSE	3053	1419	18125	72540	13	7	1065	1082	390	397	45	158	15072	71121	-	-
126	186	TRUE	FALSE	1124	386	7076	28658	12	6	1093	1092	394	397	44	158	5952	28272	-	-
127	634	TRUE	FALSE	9949	2382	46721	96214	27	8	521	506	392	388	85	156	36772	93832	-	-
128	299	TRUE	FALSE	1935	868	11503	45718	14	7	1102	1098	399	398	46	157	9568	44850	-	-
129	77	TRUE	FALSE	308	177	2233	11419	6	4	1296	1301	393	392	31	150	1925	11242	-	-
130	246	TRUE	FALSE	2272	506	14080	37406	18	6	884	881	403	406	66	156	11808	36900	-	-
131	159	TRUE	FALSE	658	374	6064	24224	9	11	1134	1137	402	395	43	161	5406	23850	-	-

132	967	TRUE	FALSE	9829	11401	39806	184494	20	34	290	302	409	396	51	213	29977	173093	-	-
133	257	TRUE	FALSE	2184	415	15291	38965	18	4	852	845	405	408	69	154	13107	38550	-	-
134	460	TRUE	FALSE	10246	1049	35086	70509	45	7	770	785	415	414	99	158	24840	69460	-	-
135	642	TRUE	FALSE	10596	1188	46548	98129	32	7	795	787	413	413	88	158	35952	96941	-	-
136	626	TRUE	FALSE	23337	1264	57141	95164	79	16	737	750	428	414	133	166	33804	93900	-	-
137	1739	TRUE	FALSE	163206	12198	255373	260872	202	13	758	773	447	425	255	156	92167	248674	-	-
138	1093	TRUE	TRUE	60683	4204	117519	164875	131	14	684	688	433	438	183	161	56836	160671	-	-
139	1107	TRUE	TRUE	77671	5169	137449	169005	173	17	710	710	434	439	227	165	59778	163836	-	-
140	743	TRUE	FALSE	26033	2211	66898	111432	76	7	662	664	431	434	131	154	40865	109221	-	-
141	293	TRUE	FALSE	1582	736	8614	42928	11	9	1239	1235	414	407	35	153	7032	42192	-	-
142	240	TRUE	FALSE	2024	523	8024	34843	17	8	1254	1253	415	406	42	151	6000	34320	-	-
143	107	TRUE	FALSE	507	309	3182	15396	11	10	1263	1264	414	407	36	151	2675	15087	-	-
144	321	TRUE	TRUE	2848	999	9268	49149	19	11	1326	1323	416	409	39	161	6420	48150	-	-
147	469	TRUE	FALSE	8117	1565	33443	69570	33	6	623	621	436	430	87	151	25326	68005	-	-
148	262	TRUE	FALSE	2867	272	13871	39572	21	2	970	965	423	416	63	152	11004	39300	-	-
149	2895	TRUE	TRUE	170968	12237	324403	434907	202	16	839	872	443	441	255	162	153435	422670	-	-
151	188	TRUE	FALSE	1484	233	9004	28433	13	4	993	987	424	424	53	154	7520	28200	-	-
153	230	TRUE	FALSE	1443	1590	8803	41380	13	12	318	309	430	420	45	185	7360	39790	-	-
155	257	TRUE	FALSE	2004	1371	10999	46089	13	43	291	297	428	421	48	217	8995	44718	-	-
156	981	TRUE	FALSE	13262	8255	43673	178949	26	22	366	344	434	422	57	196	30411	170694	-	-
157	349	TRUE	FALSE	2740	755	14606	52407	17	5	1092	1101	432	436	51	153	11866	51652	-	-
158	628	TRUE	FALSE	7040	2790	28392	108922	24	17	285	274	446	429	58	186	21352	106132	-	-
159	176	TRUE	FALSE	1465	1138	7625	31938	19	17	372	373	435	443	54	192	6160	30800	-	-
160	266	TRUE	FALSE	2690	520	12266	39888	28	6	1110	1103	434	437	64	154	9576	39368	-	-
161	319	TRUE	FALSE	3776	571	14941	47783	35	3	1127	1130	435	425	70	151	11165	47212	-	-
162	623	TRUE	FALSE	6366	3915	33155	102349	19	45	443	447	435	440	62	203	26789	98434	-	-
163	823	TRUE	FALSE	8175	8633	40272	143605	26	25	403	380	436	438	65	189	32097	134972	-	-
164	931	TRUE	FALSE	88268	5165	148783	139229	167	13	733	723	454	436	232	157	60515	134064	-	-
165	230	TRUE	FALSE	1021	405	10451	34675	9	4	1050	1047	452	451	50	153	9430	34270	-	-
167	349	TRUE	FALSE	3653	504	19707	52854	20	3	972	966	435	432	66	153	16054	52350	-	-
168	117	TRUE	FALSE	1049	560	5494	20450	17	16	277	274	446	437	55	186	4445	19890	-	-
169	153	TRUE	FALSE	2631	360	21144	23004	33	6	698	696	441	435	154	154	18513	22644	-	-
170	223	TRUE	FALSE	4574	419	21745	32977	41	8	658	661	442	443	118	154	17171	32558	-	-
171	259	TRUE	FALSE	1154	459	9960	37496	11	4	1185	1183	450	440	45	147	8806	37037	-	-
173	294	TRUE	FALSE	1454	466	9980	42214	12	8	1215	1215	451	447	41	150	8526	41748	-	-
174	509	TRUE	TRUE	8279	2743	30166	78075	55	20	1010	1011	458	459	98	168	21887	75332	-	-
175	527	TRUE	FALSE	34899	1477	68627	77365	135	10	710	716	454	445	199	154	33728	75888	-	-
176	570	TRUE	TRUE	9853	3067	35503	88567	54	25	979	979	457	459	99	175	25650	85500	-	-
177	242	TRUE	FALSE	5894	282	17510	35130	44	4	649	640	455	466	92	148	11616	34848	-	-
178	444	TRUE	FALSE	16785	1267	42093	64759	91	7	694	689	454	449	148	150	25308	63492	-	-
179	439	TRUE	TRUE	5394	1624	24710	67035	24	11	944	957	458	460	68	160	19316	65411	-	-
180	276	TRUE	FALSE	3481	311	17281	40055	24	4	619	600	454	456	74	148	13800	39744	-	-
181	328	TRUE	FALSE	6646	257	22718	47473	42	3	634	622	455	450	91	147	16072	47216	-	-
182	247	TRUE	FALSE	6751	366	19101	35934	51	3	663	656	455	450	101	147	12350	35568	-	-
183	253	TRUE	FALSE	6943	675	20605	36854	55	5	676	667	453	453	109	148	13662	36179	-	-
184	299	TRUE	FALSE	2317	550	15174	44204	17	5	1067	1063	473	469	60	151	12857	43654	-	-
185	436	TRUE	FALSE	5419	1746	23295	66710	35	8	1102	1102	475	475	76	157	17876	64964	-	-
186	531	TRUE	FALSE	8295	1773	30066	80361	51	7	1134	1129	477	476	92	155	21771	78588	-	-
188	403	TRUE	TRUE	2972	1562	16674	59594	15	9	1195	1195	479	482	49	153	13702	58032	-	-
189	386	TRUE	TRUE	2961	3561	14155	59531	13	17	1228	1230	481	488	42	162	11194	55970	-	-
190	151	TRUE	FALSE	936	670	8033	22414	11	36	607	607	478	483	58	180	7097	21744	-	-
191	492	TRUE	FALSE	7506	877	33582	70741	38	6	702	712	480	486	91	148	26076	69864	-	-
192	322	TRUE	FALSE	2167	309	18589	46354	15	2	674	666	480	473	66	145	16422	46045	-	-
193	3233	TRUE	FALSE	144219	9493	309102	465346	178	9	869	876	504	537	229	150	164883	455853	-	-
194	323	TRUE	FALSE	7056	570	25790	46759	50	8	733	742	480	476	108	151	18734	46189	-	-
199	310	TRUE	TRUE	2861	1155	17121	47035	21	8	1154	1161	500	497	67	156	14260	45880	-	-
200	915	TRUE	FALSE	18449	1212	74264	131142	48	4	781	768	505	507	109	146	55815	129930	-	-
201	1088	TRUE	FALSE	15416	2284	78520	156780	30	15	736	729	500	496	88	157	63104	154496	-	-
202	261	TRUE	TRUE	2567	1678	13007	39523	19	10	1179	1184	498	505	59	155	10440	37845	-	-
203	371	TRUE	TRUE	3057	3464	17155	57630	15	16	1192	1196	499	500	53	162	14098	54166	-	-

204	1753	TRUE	FALSE	56329	10325	131708	266263	103	29	1073	1080	511	518	146	175	75379	255938	-	-
205	353	TRUE	FALSE	3241	813	21597	53763	17	6	567	554	499	502	69	156	18356	52950	-	-
206	641	TRUE	FALSE	9998	2419	40125	99210	51	20	1108	1103	507	519	98	171	30127	96791	-	-
207	278	TRUE	TRUE	2408	1631	11026	44443	14	12	1228	1229	501	498	45	166	8618	42812	-	-
208	1045	TRUE	FALSE	8259	8231	52149	171251	14	97	454	423	516	505	56	253	43890	163020	-	-
209	995	TRUE	FALSE	10350	2064	68060	156282	17	6	497	501	529	510	75	161	57710	154218	-	-
210	727	TRUE	FALSE	14887	3529	49783	113306	55	9	1130	1131	518	529	103	160	34896	109777	-	-
211	413	TRUE	FALSE	5569	2065	26219	64015	38	10	1153	1149	511	519	88	160	20650	61950	-	-
212	400	TRUE	FALSE	3143	1751	25543	58951	17	10	695	701	517	519	73	153	22400	57200	-	-
213	408	TRUE	FALSE	2628	384	23436	61176	13	2	603	601	514	504	64	151	20808	60792	-	-
214	192	TRUE	FALSE	1031	253	11783	29053	10	2	569	565	510	508	66	152	10752	28800	-	-
215	443	TRUE	FALSE	3325	840	25918	69504	14	6	481	461	524	527	65	161	22593	68664	-	-
216	131	TRUE	FALSE	373	270	3517	19789	4	4	1330	1330	515	511	28	153	3144	19519	-	-
217	147	TRUE	TRUE	544	177	8482	22227	7	18	662	658	520	517	61	168	7938	22050	-	-
218	490	TRUE	FALSE	5834	2140	32784	72210	26	10	764	760	528	535	81	153	26950	70070	-	-
219	349	TRUE	FALSE	4410	874	20813	54620	25	6	1175	1172	537	545	72	160	16403	53746	-	-
220	563	TRUE	FALSE	13607	2099	45135	82045	70	10	806	814	527	541	126	152	31528	79946	-	-
221	354	TRUE	FALSE	4732	496	26326	53950	26	4	542	529	532	526	87	155	21594	53454	-	-
222	416	TRUE	FALSE	4806	478	28518	62878	32	5	557	566	534	529	89	155	23712	62400	-	-
223	618	TRUE	FALSE	4750	973	35032	93673	15	5	955	967	532	537	64	155	30282	92700	-	-
224	620	TRUE	FALSE	11307	2057	42307	96297	61	13	1016	1008	535	548	111	165	31000	94240	-	-
225	504	TRUE	TRUE	8373	1892	35085	80012	38	7	1152	1163	543	545	91	162	26712	78120	-	-
226	646	TRUE	TRUE	18454	1707	48816	101191	96	7	1077	1084	536	526	143	161	30362	99484	-	-
227	399	TRUE	FALSE	3644	378	22796	61026	24	3	987	996	539	545	72	155	19152	60648	-	-
229	309	TRUE	FALSE	3610	264	19369	48133	27	3	1122	1121	545	537	78	158	15759	47869	-	-
230	345	TRUE	FALSE	4296	445	23961	54265	29	4	1132	1138	540	548	86	160	19665	53820	-	-
232	730	TRUE	FALSE	12799	1815	54409	110585	48	7	719	716	550	542	105	156	41610	108770	-	-
233	729	TRUE	FALSE	17886	4234	63813	113584	73	21	532	528	551	554	136	171	45927	109350	-	-
235	652	TRUE	TRUE	10708	2254	57000	102010	45	9	518	515	552	553	116	162	46292	99756	-	-
236	532	TRUE	FALSE	4792	733	36180	80001	21	15	689	683	551	543	80	164	31388	79268	-	-
239	565	TRUE	FALSE	7179	702	18479	84887	30	3	1338	1341	565	552	50	152	11300	84185	-	-
240	478	TRUE	TRUE	7239	2904	31617	76994	32	15	1119	1118	555	562	83	170	24378	74090	-	-
242	328	TRUE	FALSE	1886	1301	17630	52141	17	9	1046	1056	559	560	65	164	15744	50840	-	-
243	245	TRUE	FALSE	2220	754	14470	38728	17	6	1078	1078	551	561	67	161	12250	37974	-	-
244	327	TRUE	FALSE	3965	456	20969	49506	23	4	570	565	561	556	75	154	17004	49050	-	-
245	569	TRUE	FALSE	9290	733	42292	86083	56	3	816	804	564	563	114	153	33002	85350	-	-
246	280	TRUE	TRUE	2767	2232	15647	45632	18	23	1098	1101	561	564	64	178	12880	43400	-	-
247	625	TRUE	FALSE	7872	542	45372	94917	43	3	780	779	565	558	103	154	37500	94375	-	-
249	229	TRUE	FALSE	2352	252	16550	34831	22	11	745	745	566	557	84	162	14198	34579	-	-
250	2301	TRUE	FALSE	65432	5036	191987	357089	77	7	495	488	588	601	132	160	126555	352053	-	-
255	1334	TRUE	FALSE	38158	2262	111528	203696	89	5	727	727	585	595	144	156	73370	201434	-	-
256	369	TRUE	FALSE	3276	327	22464	57105	24	6	1066	1072	573	581	76	160	19188	56778	-	-
257	669	TRUE	FALSE	6509	1340	42635	101690	17	5	878	877	580	592	71	155	36126	100350	-	-
258	276	TRUE	FALSE	2161	605	17065	42281	13	5	913	917	579	580	67	156	14904	41676	-	-
259	562	TRUE	FALSE	7214	1607	38123	88717	44	8	1123	1137	577	586	99	163	30909	87110	-	-
260	514	TRUE	FALSE	4367	1028	33151	78108	16	7	602	612	589	577	72	157	28784	77080	-	-
262	196	TRUE	FALSE	1237	189	10841	30373	14	2	1083	1080	581	580	63	156	9604	30184	-	-
263	44	TRUE	FALSE	70	87	1830	7347	3	4	280	278	576	574	43	169	1760	7260	-	-
264	487	TRUE	FALSE	9216	770	35514	74307	52	6	576	575	588	600	106	157	26298	73537	-	-
265	200	TRUE	FALSE	569	377	7969	32177	5	3	362	362	584	577	42	162	7400	31800	-	-
267	378	TRUE	TRUE	4116	4339	19614	66709	18	27	1226	1231	589	592	59	192	15498	62370	-	-
268	372	TRUE	FALSE	3109	953	17245	61217	17	14	274	278	597	582	55	176	14136	60264	-	-
269	139	TRUE	FALSE	635	382	8419	23595	10	5	1159	1161	586	583	66	172	7784	23213	-	-
271	210	TRUE	FALSE	1239	527	14049	32237	12	5	843	843	591	596	73	156	12810	31710	-	-
272	2284	TRUE	TRUE	157330	9274	278382	361009	202	37	992	979	606	624	255	191	121052	351735	-	-
273	1728	TRUE	FALSE	87028	4346	182068	266996	188	10	919	920	615	612	243	162	95040	262650	-	-
274	462	TRUE	FALSE	4576	1140	29524	72750	19	5	1109	1108	603	591	73	160	24948	71610	-	-
275	333	TRUE	FALSE	2112	1202	14433	53816	12	8	1258	1254	610	603	49	166	12321	52614	-	-
276	1141	TRUE	FALSE	22836	2499	84450	178213	67	9	1076	1074	614	601	121	163	61614	175714	-	-
277	291	TRUE	FALSE	2609	1203	19487	48927	17	13	1147	1145	604	596	75	177	16878	47724	-	-

278	555	TRUE	FALSE	7652	1631	37622	85991	30	6	948	946	613	622	84	158	29970	84360	-	-
279	609	TRUE	FALSE	6068	3167	38345	104261	22	11	1183	1183	607	606	75	177	32277	101094	-	-
280	1662	TRUE	FALSE	43046	4591	136118	262201	76	25	489	477	627	629	132	180	93072	257610	-	-
281	967	TRUE	FALSE	30049	1310	79366	148287	86	5	849	840	617	625	137	157	49317	146977	-	-
282	857	TRUE	FALSE	18880	905	63444	131158	49	5	874	875	613	626	101	157	44564	130253	-	-
283	316	TRUE	FALSE	2845	1567	21489	51179	21	12	1128	1139	609	610	80	169	18644	49612	-	-
284	579	TRUE	FALSE	15376	500	44905	88502	65	4	821	818	615	602	116	156	29529	88002	-	-
285	584	TRUE	FALSE	15144	684	44928	89451	69	13	808	801	616	605	120	165	29784	88767	-	-
286	196	TRUE	FALSE	1717	458	13085	33190	15	5	1165	1158	610	611	73	172	11368	32732	-	-
288	319	TRUE	FALSE	3600	294	20826	49420	26	2	1021	1021	615	615	80	156	17226	49126	-	-
289	1029	TRUE	FALSE	13949	1423	65399	157831	28	4	770	736	617	619	78	156	51450	156408	-	-
290	884	TRUE	TRUE	13576	8088	47168	150412	21	28	1210	1227	627	623	59	189	33592	142324	-	-
292	551	TRUE	FALSE	4475	1232	40290	86034	17	6	560	554	615	624	82	160	35815	84802	-	-
293	411	TRUE	FALSE	2654	620	26903	62681	13	4	670	669	622	620	72	155	24249	62061	-	-
294	728	TRUE	FALSE	8532	3898	34740	118194	30	12	1251	1238	631	620	66	169	26208	114296	-	-
296	372	TRUE	TRUE	4697	2096	24785	60872	31	16	1130	1132	629	629	85	174	20088	58776	-	-
297	747	TRUE	TRUE	14618	6425	49727	127439	39	24	1182	1183	628	635	86	186	35109	121014	-	-
298	303	TRUE	FALSE	3240	1057	19299	50143	20	7	1168	1168	626	618	73	169	16059	49086	-	-
299	298	TRUE	FALSE	2287	602	21359	45600	16	5	627	620	627	635	80	156	19072	44998	-	-
301	834	TRUE	FALSE	37786	2004	84490	132942	157	6	1083	1085	646	646	213	163	46704	130938	-	-
302	563	TRUE	FALSE	6330	1964	41236	88666	23	8	542	533	642	635	85	162	34906	86702	-	-
303	137	TRUE	FALSE	398	211	5878	21720	5	3	369	366	627	623	45	160	5480	21509	-	-
304	412	TRUE	FALSE	9255	2638	31503	67322	47	34	471	468	638	643	101	191	22248	64684	-	-
305	1424	TRUE	FALSE	49845	5159	129589	224455	116	8	1028	1046	652	636	172	162	79744	219296	-	-
306	368	TRUE	FALSE	2973	945	25789	57617	15	17	526	521	639	648	77	171	22816	56672	-	-
307	985	TRUE	FALSE	16942	2898	71117	154588	40	17	969	959	652	650	95	171	54175	151690	-	-
308	983	TRUE	TRUE	21035	4550	80015	153966	73	28	647	647	648	651	133	180	58980	149416	-	-
309	562	TRUE	TRUE	13931	4070	43717	93428	74	19	1132	1133	641	647	127	178	29786	89358	-	-
310	822	TRUE	TRUE	12791	3628	62933	130216	50	19	607	606	648	650	111	173	50142	126588	-	-
311	244	TRUE	FALSE	2026	715	11298	39511	18	12	370	370	645	648	56	171	9272	38796	-	-
312	198	TRUE	FALSE	1889	418	13769	31306	17	5	565	562	644	644	77	161	11880	30888	-	-
313	496	TRUE	FALSE	5610	1527	34378	77911	25	8	571	572	646	663	83	162	28768	76384	-	-
314	622	TRUE	TRUE	36312	1785	74254	100061	162	39	1084	1094	655	651	223	197	37942	98276	-	-
315	493	TRUE	FALSE	3553	1682	32640	77111	14	14	948	944	648	646	73	167	29087	75429	-	-
316	256	TRUE	FALSE	1494	256	15062	39424	11	3	833	832	649	657	64	156	13568	39168	-	-
317	945	TRUE	TRUE	26502	3946	78477	153256	83	14	1109	1100	671	670	138	172	51975	149310	-	-
318	535	TRUE	TRUE	10546	15782	38366	99242	55	91	469	467	651	651	107	247	27820	83460	-	-
319	1104	TRUE	TRUE	18197	8072	81125	175880	42	27	708	704	673	675	99	179	62928	167808	-	-
320	502	TRUE	TRUE	6597	6131	24167	85447	26	29	370	370	655	655	61	187	17570	79316	-	-
321	610	TRUE	FALSE	4339	4907	36669	98237	13	29	890	884	657	670	66	182	32330	93330	-	-
322	214	TRUE	TRUE	4049	656	16247	34682	36	7	1126	1131	653	656	93	166	12198	34026	-	-
323	359	TRUE	FALSE	3789	1175	26047	56820	21	13	997	1007	661	664	83	168	22258	55645	-	-
324	550	TRUE	FALSE	4261	1144	15261	85844	12	4	1331	1337	658	650	32	158	11000	84700	-	-
325	198	TRUE	TRUE	3376	352	14464	31636	36	7	1136	1133	656	657	92	165	11088	31284	-	-
326	844	TRUE	TRUE	42589	5772	94073	140812	184	26	1070	1072	671	668	245	186	51484	135040	-	-
327	869	TRUE	TRUE	12366	8633	64506	140721	33	26	663	674	674	676	93	178	52140	132088	-	-
328	272	TRUE	FALSE	3938	1367	18898	43255	30	11	983	986	667	670	85	165	14960	41888	-	-
330	285	TRUE	FALSE	4947	1454	19767	45627	39	16	972	962	661	671	91	171	14820	44173	-	-
331	356	TRUE	TRUE	3685	1743	23621	58347	19	18	1172	1180	671	675	75	177	19936	56604	-	-
332	1804	TRUE	TRUE	48666	20027	122630	299647	62	32	419	392	682	690	103	187	73964	279620	-	-
333	1240	TRUE	TRUE	35864	6778	106544	200218	116	13	1028	1030	670	674	173	169	70680	193440	-	-
334	406	TRUE	FALSE	5186	1185	23456	64927	25	6	1209	1206	673	674	70	163	18270	63742	-	-
335	611	TRUE	FALSE	10328	1073	29880	96388	39	4	1259	1273	676	669	71	160	19552	95315	-	-
336	479	TRUE	TRUE	7202	4794	31631	79039	31	25	440	443	680	686	82	180	24429	74245	-	-
337	401	TRUE	TRUE	5260	2486	29320	65844	23	12	1137	1132	672	674	83	170	24060	63358	-	-
338	1693	TRUE	TRUE	113439	9757	209940	270479	198	24	593	606	693	675	255	178	96501	260722	-	-
339	260	TRUE	FALSE	2432	902	17772	40682	22	10	997	992	678	673	81	163	15340	39780	-	-
340	794	TRUE	FALSE	26418	4708	78822	126190	108	28	620	620	693	676	174	181	52404	121482	-	-
341	504	TRUE	FALSE	3865	1257	29065	78369	15	13	879	879	687	675	65	166	25200	77112	-	-
342	816	TRUE	TRUE	29019	3623	77163	129287	111	19	570	576	695	679	170	173	48144	125664	-	-

343	1135	TRUE	TRUE	58377	6558	124207	185888	137	17	1187	1183	713	676	195	175	65830	179330	-	-
344	1873	TRUE	TRUE	59735	6527	164623	304334	86	9	1121	1130	710	679	142	168	104888	297807	-	-
345	940	TRUE	TRUE	34983	2696	91383	153089	87	9	1167	1164	701	706	147	169	56400	150393	-	-
346	1078	TRUE	FALSE	50156	2248	101900	171491	106	9	1209	1223	703	722	154	166	51744	169243	-	-
347	673	TRUE	FALSE	6558	652	39535	103616	18	13	914	900	692	687	67	166	32977	102964	-	-
348	552	TRUE	FALSE	5362	1276	40138	86836	25	16	552	546	693	681	88	171	34776	85560	-	-
349	916	TRUE	FALSE	17666	2639	68962	147367	78	6	1071	1065	684	683	134	164	51296	144728	-	-
350	2074	TRUE	FALSE	110747	6487	177115	325879	165	12	1260	1239	716	714	197	166	66368	319392	-	-
351	1090	TRUE	FALSE	13182	2762	67682	171709	22	18	498	464	698	685	72	173	54500	168947	-	-
352	286	TRUE	FALSE	1661	265	16819	44023	12	2	983	981	695	690	65	155	15158	43758	-	-
353	557	TRUE	FALSE	7811	1225	37332	88674	26	7	1054	1063	706	694	79	164	29521	87449	-	-
354	194	TRUE	FALSE	1365	247	13199	30317	13	6	505	504	698	689	74	161	11834	30070	-	-
355	253	TRUE	FALSE	1383	449	14539	39157	10	10	950	951	697	691	62	163	13156	38708	-	-
356	189	TRUE	FALSE	1454	223	13172	29518	13	3	521	512	699	691	75	158	11718	29295	-	-
357	553	TRUE	TRUE	5641	6226	37715	92494	19	29	728	730	708	711	77	185	32074	86268	-	-
358	271	TRUE	FALSE	863	548	15497	42011	5	5	800	806	699	698	59	158	14634	41463	-	-
359	470	TRUE	FALSE	6108	690	25378	74010	22	10	398	389	708	696	63	166	19270	73320	-	-
360	716	TRUE	TRUE	13313	2829	57705	113093	48	8	648	655	714	710	110	162	44392	110264	-	-
361	881	TRUE	FALSE	18564	1696	47637	136489	48	4	1281	1296	717	708	81	157	29073	134793	-	-
362	305	TRUE	FALSE	3053	412	10983	47077	22	3	1327	1334	710	708	48	156	7930	46665	-	-
363	527	TRUE	TRUE	4759	4615	35325	85773	16	20	689	691	711	711	74	174	30566	81158	-	-
364	171	TRUE	FALSE	3563	205	15875	27907	32	3	1152	1155	711	717	104	165	12312	27702	-	-
365	111	TRUE	FALSE	3178	397	10837	17824	48	8	1225	1225	715	721	117	165	7659	17427	-	-
366	875	TRUE	FALSE	40855	2460	97730	137202	140	6	618	630	717	719	205	160	56875	134742	-	-
367	428	TRUE	FALSE	4422	456	24966	65512	24	3	930	928	714	705	72	155	20544	65056	-	-
368	406	TRUE	FALSE	4834	437	17420	62555	23	3	1323	1317	716	711	54	156	12586	62118	-	-
369	487	TRUE	FALSE	5756	1270	23775	77242	23	9	361	369	719	730	60	165	18019	75972	-	-
370	759	TRUE	FALSE	6515	2619	57368	120264	20	8	513	507	729	727	87	163	50853	117645	-	-
371	555	TRUE	FALSE	6508	772	32593	85129	23	2	973	975	718	707	70	154	26085	84357	-	-
372	39	TRUE	FALSE	380	48	1784	6015	15	2	1313	1314	718	709	51	155	1404	5967	-	-
373	981	TRUE	FALSE	21269	2907	70319	154962	51	8	1063	1062	725	733	101	163	49050	152055	-	-
374	107	TRUE	FALSE	3857	321	12203	17441	63	6	1204	1202	715	725	141	166	8346	17120	-	-
375	430	TRUE	FALSE	5944	1117	26154	68627	24	7	408	411	722	727	71	164	20210	67510	-	-
376	129	TRUE	FALSE	506	109	7085	20100	8	3	844	851	718	719	59	158	6579	19991	-	-
377	537	TRUE	FALSE	4316	2208	38684	85443	14	8	493	475	723	726	78	163	34368	83235	-	-
378	568	TRUE	TRUE	8401	1242	33961	88710	31	13	1015	1021	724	735	76	167	25560	87468	-	-
379	207	TRUE	TRUE	1180	731	13186	32402	11	8	726	722	721	717	69	161	12006	31671	-	-
380	691	TRUE	FALSE	4177	1348	38036	107071	13	8	843	848	726	725	62	161	33859	105723	-	-
381	410	TRUE	TRUE	10630	1240	35230	66840	69	9	1153	1162	727	729	129	169	24600	65600	-	-
382	912	TRUE	TRUE	23639	4215	71975	147399	73	12	1114	1117	733	746	126	169	48336	143184	-	-
383	347	TRUE	FALSE	7302	442	28816	55615	50	3	1132	1133	732	723	112	162	21514	55173	-	-
384	609	TRUE	FALSE	9551	2121	43046	99561	30	9	421	438	739	731	85	169	33495	97440	-	-
386	1030	TRUE	FALSE	35753	1426	99612	159016	105	4	642	644	752	725	105	157	63859	157590	-	-
387	1022	TRUE	TRUE	28019	12832	80141	176352	82	45	1220	1224	734	747	133	205	52122	163520	-	-
388	1044	TRUE	FALSE	17097	1101	41109	160823	42	4	1328	1322	739	742	65	157	24012	159722	-	-
389	695	TRUE	TRUE	14223	7536	54533	116651	52	43	454	452	746	749	110	200	40310	109115	-	-
390	81	TRUE	FALSE	406	100	4213	12493	10	3	926	927	727	734	57	156	3807	12393	-	-
391	544	TRUE	TRUE	15206	1897	48390	91113	60	12	1192	1181	735	746	121	176	33184	89216	-	-
392	1024	TRUE	FALSE	40619	5107	74411	163824	93	17	1268	1240	738	745	126	172	33792	158717	-	-
393	4152	TRUE	FALSE	250785	14231	387801	666095	222	15	357	370	771	775	255	172	137016	651864	-	-
394	289	TRUE	FALSE	3659	1044	18398	46995	27	8	413	407	736	732	78	167	14739	45951	-	-
395	726	TRUE	FALSE	15297	1574	37803	113378	65	13	1283	1295	739	738	96	167	22506	111804	-	-
396	506	TRUE	FALSE	11153	854	44549	78272	51	10	617	604	752	741	117	163	33396	77418	-	-
397	343	TRUE	FALSE	1840	393	17961	52186	11	4	918	914	742	739	58	155	16121	51793	-	-
398	918	TRUE	FALSE	11060	1824	71648	144114	30	9	483	483	747	745	96	164	60588	142290	-	-
399	319	TRUE	FALSE	2505	697	18774	49185	18	4	878	876	743	735	69	156	16269	48488	-	-
400	648	TRUE	FALSE	11169	752	50697	99896	37	3	673	672	753	738	98	156	39528	99144	-	-
401	426	TRUE	FALSE	4376	682	26528	65433	31	3	819	819	750	737	83	155	22152	64751	-	-
402	679	TRUE	FALSE	6059	776	50194	105340	21	5	591	594	751	755	86	159	44135	104564	-	-
403	346	TRUE	FALSE	2169	506	20853	53444	13	5	778	771	751	749	67	158	18684	52938	-	-

404	542	TRUE	TRUE	7862	1175	34420	84641	41	9	1028	1027	755	760	90	163	26558	83466	-	-
405	30	TRUE	FALSE	84	106	2274	4966	5	6	430	433	743	747	78	168	2190	4860	-	-
406	928	TRUE	TRUE	22167	21796	74135	167492	67	98	451	454	757	755	123	255	51968	145696	-	-
407	454	TRUE	FALSE	4733	1484	28341	70492	23	78	765	753	762	761	75	230	23608	69008	-	-
408	1017	TRUE	FALSE	12563	2946	70532	156513	25	5	705	708	772	761	82	156	57969	153567	-	-
409	469	TRUE	FALSE	3357	3075	18365	77177	15	18	274	274	761	767	47	176	15008	74102	-	-
410	1262	TRUE	TRUE	28811	4193	93173	197279	71	62	1077	1077	780	788	122	215	64362	193086	-	-
411	270	TRUE	FALSE	2192	1390	10292	43780	13	11	291	292	762	769	43	168	8100	42390	-	-
412	536	TRUE	FALSE	10170	1315	45010	83323	50	6	616	617	760	762	115	159	34840	82008	-	-
413	1037	TRUE	TRUE	36572	3545	99829	161169	106	12	643	643	760	765	167	164	63257	157624	-	-
414	344	TRUE	FALSE	3360	405	19184	52349	21	3	930	935	769	759	67	154	15824	51944	-	-
415	1236	TRUE	TRUE	33633	3850	107793	192958	99	9	1126	1130	782	783	159	162	74160	189108	-	-
416	544	TRUE	TRUE	6175	2871	38814	86647	25	15	537	539	777	780	85	169	32639	83776	-	-
419	936	TRUE	FALSE	42221	2324	96509	151147	111	5	1217	1219	812	765	169	164	54288	148823	-	-
420	1357	TRUE	FALSE	58455	5314	137161	218363	115	23	1202	1199	810	767	173	180	78706	213049	-	-
421	574	TRUE	FALSE	13166	974	47606	88222	52	5	670	662	776	777	112	157	34440	87248	-	-
422	1273	TRUE	FALSE	32051	6991	94428	206852	77	10	382	409	789	779	126	167	62377	199861	-	-
423	500	TRUE	FALSE	3863	1042	32863	78042	16	5	564	565	778	780	74	159	29000	77000	-	-
424	840	TRUE	FALSE	17659	3047	69739	133247	59	7	1170	1169	792	773	121	162	52080	130200	-	-
425	1456	TRUE	FALSE	22960	3911	54992	223767	38	5	1338	1332	814	804	60	156	32032	219856	-	-
427	498	TRUE	FALSE	5421	1830	34803	79518	17	9	460	456	781	776	76	165	29382	77688	-	-
428	221	TRUE	FALSE	962	407	12675	33557	7	8	827	819	774	777	60	158	11713	33150	-	-
429	265	TRUE	FALSE	1555	300	14805	40315	14	3	973	971	780	780	64	154	13250	40015	-	-
430	462	TRUE	FALSE	4236	915	26412	69753	19	6	907	906	789	775	67	155	22176	68838	-	-
431	338	TRUE	FALSE	2670	302	21598	52016	14	2	592	593	782	775	70	155	18928	51714	-	-
432	314	TRUE	FALSE	1371	343	18327	48071	7	5	747	745	778	782	61	157	16956	47728	-	-
433	292	TRUE	FALSE	1675	744	15399	44252	11	4	946	939	783	778	58	153	13724	43508	-	-
435	407	TRUE	FALSE	4615	590	26186	62454	26	3	635	631	795	794	79	155	21571	61864	-	-
436	677	TRUE	FALSE	15463	608	49990	103512	72	14	667	684	795	803	123	166	34527	102904	-	-
437	320	TRUE	FALSE	3584	270	20864	48910	21	3	706	707	794	790	75	155	17280	48640	-	-
438	700	TRUE	FALSE	21314	1020	60514	107420	63	4	1086	1097	806	790	119	156	39200	106400	-	-
439	1271	TRUE	FALSE	31887	3554	86540	192933	65	13	989	978	822	821	108	162	54653	189379	-	-
440	1244	TRUE	FALSE	25073	2449	58661	194025	49	34	1286	1299	814	830	76	188	33588	191576	-	-
441	416	TRUE	FALSE	9896	337	37768	63985	45	2	1107	1103	812	812	112	155	27872	63648	-	-
442	1122	TRUE	FALSE	36851	2368	121001	176277	82	6	1171	1182	809	798	157	161	84150	173909	-	-
443	1142	TRUE	FALSE	30847	2917	110787	176501	70	8	1130	1151	809	799	140	160	79940	173584	-	-
444	1134	TRUE	FALSE	54746	2848	120517	175216	139	10	531	536	820	815	197	162	65771	172368	-	-
445	1840	TRUE	FALSE	41391	6987	144431	294027	56	19	429	402	820	824	112	175	103040	287040	-	-
446	510	TRUE	FALSE	14640	8263	47280	78642	58	20	511	511	820	805	122	158	32640	70379	-	-
447	555	TRUE	FALSE	9113	563	37973	85478	44	4	552	559	818	804	96	157	28860	84915	-	-
448	481	TRUE	FALSE	10538	974	35550	73124	61	4	1000	1010	821	815	113	154	25012	72150	-	-
449	1740	TRUE	FALSE	57237	23452	159897	268792	92	21	482	498	827	824	151	162	102660	245340	-	-
450	290	TRUE	FALSE	2854	181	17644	44249	20	4	716	716	811	805	71	156	14790	44068	-	-
453	445	TRUE	FALSE	4194	737	27779	68377	22	6	619	620	816	821	75	158	23585	67640	-	-
454	293	TRUE	FALSE	1594	304	17123	45133	10	3	590	589	814	808	63	156	15529	44829	-	-
455	293	TRUE	FALSE	2527	475	16298	43546	17	4	937	942	821	820	64	151	13771	43071	-	-
456	307	TRUE	FALSE	4145	561	24407	48146	29	14	464	465	827	827	95	169	20262	47585	-	-
458	763	TRUE	FALSE	21135	1521	74545	116734	58	4	1100	1110	829	822	128	155	53410	115213	-	-
459	1229	TRUE	FALSE	49895	3955	127322	189534	96	6	1142	1138	847	848	159	157	77427	185579	-	-
461	386	TRUE	FALSE	14075	621	45341	61223	76	4	1202	1202	843	828	157	161	31266	60602	-	-
462	666	TRUE	FALSE	6709	910	40009	102142	26	3	612	604	845	838	76	155	33300	101232	-	-
463	557	TRUE	FALSE	32799	965	63434	88414	149	3	1233	1224	843	829	204	160	30635	87449	-	-
465	492	TRUE	FALSE	3502	752	26134	75040	14	3	751	755	840	834	60	154	22632	74288	-	-
466	201	TRUE	FALSE	1294	274	11545	30826	17	3	574	579	846	836	64	155	10251	30552	-	-
467	1567	TRUE	FALSE	32505	2000	103020	233913	63	3	881	860	857	846	112	151	70515	231913	-	-
468	725	TRUE	FALSE	32457	1529	70157	114629	166	8	1233	1243	851	854	218	164	37700	113100	-	-
469	607	TRUE	FALSE	14719	1024	46283	93288	58	4	536	526	837	837	110	156	31564	92264	-	-
470	1840	TRUE	FALSE	41756	8983	130076	294183	44	13	427	432	852	860	12	168	88320	285200	-	-
471	702	TRUE	FALSE	15174	1550	57294	106850	57	3	1101	1121	848	838	117	153	42120	105300	-	-
472	620	TRUE	FALSE	6390	454	37390	94690	22	2	694	696	852	843	72	154	31000	94236	-	-

474	311	TRUE	FALSE	2571	263	17188	47223	16	7	779	778	857	847	63	158	14617	46960	-	-
475	83	TRUE	FALSE	148	65	4630	12681	3	2	654	657	847	853	57	154	4482	12616	-	-
476	779	TRUE	FALSE	10707	1081	36414	121826	22	4	1287	1288	872	853	55	159	25707	120745	-	-
477	648	TRUE	FALSE	8590	503	39046	98996	35	3	579	571	863	854	82	155	30456	98493	-	-
478	120	TRUE	FALSE	2165	293	15845	18893	34	3	1202	1215	850	849	148	158	13680	18600	-	-
480	589	TRUE	FALSE	8745	1074	34661	88835	36	13	997	994	865	855	80	162	25916	87761	-	-
481	397	TRUE	FALSE	3399	632	19279	58991	14	4	917	917	864	867	54	151	15880	58359	-	-
482	303	TRUE	FALSE	11988	482	29865	47144	80	4	1204	1205	858	857	139	158	17877	46662	-	-
483	475	TRUE	FALSE	30019	1193	62319	74818	145	3	1233	1220	861	857	213	158	32300	73625	-	-
485	1924	TRUE	FALSE	86526	5712	155790	300084	158	8	1234	1268	873	872	194	161	69264	294372	-	-
486	181	TRUE	FALSE	5765	393	16806	27905	70	4	1148	1146	860	860	131	156	11041	27512	-	-
487	320	TRUE	FALSE	1901	470	11180	52310	12	3	312	313	871	868	41	165	9279	51840	-	-
488	441	TRUE	FALSE	3417	767	25467	67358	16	3	704	711	873	867	66	154	22050	66591	-	-
490	394	TRUE	FALSE	4653	441	21595	59935	28	3	748	739	873	875	71	154	16942	59494	-	-
491	632	TRUE	FALSE	7526	1669	27118	104053	23	7	327	331	876	881	54	169	19592	102384	-	-
492	519	TRUE	FALSE	15523	1664	45106	80033	74	6	1148	1158	869	870	131	157	29583	78369	-	-
493	133	TRUE	FALSE	755	121	3814	20470	8	2	1339	1338	868	869	31	155	3059	20349	-	-
494	794	TRUE	FALSE	6627	1018	52679	122499	19	12	540	547	899	887	77	165	46052	121481	-	-
495	125	TRUE	FALSE	986	121	4736	19496	16	11	1298	1301	873	882	46	166	3750	19375	-	-
496	412	TRUE	FALSE	4672	1481	22800	62869	27	9	996	989	880	878	71	158	18128	61388	-	-
497	878	TRUE	FALSE	13084	3816	59618	139906	28	22	473	470	889	897	81	177	46534	136090	-	-
498	1223	TRUE	FALSE	26609	7411	77975	203091	59	16	403	399	896	898	101	176	51366	195680	-	-
499	335	TRUE	FALSE	3728	630	21818	50880	27	7	1098	1100	882	888	81	157	18090	50250	-	-
500	467	TRUE	FALSE	9607	1503	36226	72487	65	12	1149	1160	885	892	122	164	26619	70984	-	-
501	947	TRUE	FALSE	7495	2475	52004	146419	16	6	705	719	897	903	63	158	44509	143944	-	-
502	187	TRUE	FALSE	761	305	6371	30973	8	4	288	285	893	894	38	168	5610	30668	-	-
503	398	TRUE	FALSE	2195	932	19707	61428	16	4	739	735	900	891	60	156	17512	60496	-	-
504	2213	TRUE	FALSE	17766	7143	64239	341306	16	28	1339	1348	969	944	37	179	46473	334163	-	-
505	477	TRUE	FALSE	5253	579	28626	73560	18	3	561	565	912	910	67	156	23373	72981	-	-
506	860	TRUE	FALSE	8859	2335	49279	129615	19	6	967	976	917	919	66	154	40420	127280	-	-
507	692	TRUE	FALSE	8882	1563	42098	104669	28	6	1051	1036	921	922	76	155	33216	103106	-	-
508	1069	TRUE	FALSE	14925	1441	77996	163928	42	12	1124	1120	924	934	101	164	63071	162487	-	-
509	1024	TRUE	FALSE	14031	2200	60111	153752	36	8	899	891	924	914	81	156	46080	151552	-	-
510	164	TRUE	FALSE	1184	296	11024	25388	16	3	540	544	916	912	76	156	9840	25092	-	-
511	686	TRUE	FALSE	9607	1103	47337	104003	29	5	1088	1085	921	926	84	155	37730	102900	-	-
512	628	TRUE	FALSE	4990	1539	35762	94483	16	5	852	850	925	935	65	153	30772	92944	-	-
513	980	TRUE	FALSE	11377	2717	53517	159517	22	9	403	387	926	912	65	169	42140	156800	-	-
514	1583	TRUE	FALSE	49807	3026	127374	249974	77	6	466	450	940	935	126	162	77567	246948	-	-
515	472	TRUE	FALSE	7155	873	30755	73558	40	5	613	624	927	920	90	159	23600	72685	-	-
516	2989	TRUE	FALSE	35674	4209	247893	458502	36	6	1155	1157	980	948	107	158	212219	454293	-	-
517	248	TRUE	FALSE	2065	556	16697	38500	16	4	540	546	919	925	75	157	14632	37944	-	-
518	203	TRUE	FALSE	593	154	10336	30806	6	2	788	781	924	923	54	153	9743	30652	-	-
519	1297	TRUE	FALSE	16909	2534	97323	199678	33	4	1214	1192	951	922	95	156	80414	197144	-	-
520	157	TRUE	FALSE	631	289	7696	23839	9	4	1017	1023	931	935	54	154	7065	23550	-	-
521	86	TRUE	FALSE	164	74	3518	14608	3	3	342	343	933	931	42	172	3354	14534	-	-
522	1152	TRUE	TRUE	41493	4805	100245	181032	145	20	610	611	950	952	196	173	58752	176227	-	-
523	901	TRUE	FALSE	12502	1685	55750	136835	26	6	962	957	945	938	74	156	43248	135150	-	-
524	156	TRUE	FALSE	2085	129	9573	23685	22	2	1052	1056	935	932	70	153	7488	23556	-	-
525	908	TRUE	FALSE	23369	2294	68769	142126	77	9	574	578	951	953	127	163	45400	139832	-	-
526	750	TRUE	FALSE	12350	869	55100	114119	39	13	1084	1095	947	937	96	164	42750	113250	-	-
527	1199	TRUE	TRUE	12897	2357	99225	184604	30	9	1196	1206	974	971	102	161	86328	182247	-	-
528	230	TRUE	FALSE	1782	373	17192	35333	18	3	1103	1114	950	948	85	155	15410	34960	-	-
529	587	TRUE	TRUE	7827	1388	37177	91199	32	10	644	637	950	950	82	163	29350	89811	-	-
530	778	TRUE	FALSE	13099	1768	48109	120024	40	5	749	760	952	942	85	157	35010	118256	-	-
531	628	TRUE	FALSE	10705	853	40221	95681	52	3	994	1006	957	951	99	154	29516	94828	-	-
532	522	TRUE	FALSE	7326	1094	33948	81482	27	4	548	547	951	954	78	158	26622	80388	-	-
533	762	TRUE	FALSE	7354	917	46216	117497	17	11	700	690	955	958	68	164	38862	116580	-	-
534	305	TRUE	FALSE	5000	400	19945	46455	32	5	1053	1053	949	952	81	156	14945	46055	-	-
535	333	TRUE	FALSE	2247	513	19563	53460	12	6	418	405	954	945	64	165	17316	52947	-	-
536	3147	TRUE	TRUE	174698	7792	338342	498724	162	8	476	458	974	984	214	164	163644	490932	-	-

537	1337	TRUE	TRUE	32730	4495	95569	206382	77	8	977	963	983	971	124	159	62839	201887	-	-
539	381	TRUE	FALSE	2928	686	21216	57836	16	5	922	920	960	965	64	155	18288	57150	-	-
540	466	TRUE	FALSE	6462	2031	37218	72863	26	21	1105	1108	959	967	92	173	30756	70832	-	-
541	486	TRUE	FALSE	4286	520	27614	73906	16	3	876	873	963	955	64	154	23328	73386	-	-
542	268	TRUE	FALSE	1974	394	14838	40862	12	4	798	798	961	960	60	155	12864	40468	-	-
543	543	TRUE	FALSE	10414	420	37021	82952	43	3	1052	1049	962	970	92	155	26607	82532	-	-
545	927	TRUE	FALSE	5517	3791	39816	154892	11	7	334	330	973	981	48	170	34299	151101	-	-
547	201	TRUE	FALSE	3264	1457	16731	32210	25	24	1214	1216	963	969	92	177	13467	30753	-	-
548	496	TRUE	FALSE	5962	894	31258	75789	31	5	823	819	975	967	82	156	25296	74895	-	-
549	851	TRUE	TRUE	12076	4946	66540	135149	32	32	1104	1108	975	971	96	185	54464	130203	-	-
551	1976	TRUE	FALSE	53799	3485	154575	301861	95	5	786	807	1008	1013	146	156	100776	298376	-	-
552	403	TRUE	FALSE	3215	794	22156	62044	14	7	1250	1255	978	971	61	159	18941	61250	-	-
553	830	TRUE	TRUE	16403	2367	62883	126867	52	8	674	682	994	976	108	158	46480	124500	-	-
556	501	TRUE	FALSE	6519	756	32070	76908	21	11	1059	1052	974	991	72	163	25551	76152	-	-
557	192	TRUE	FALSE	2773	690	14485	29682	34	6	699	693	989	977	95	157	11712	28992	-	-
558	883	TRUE	FALSE	7835	1762	61698	135095	17	4	634	632	995	981	78	155	53863	133333	-	-
559	227	TRUE	FALSE	2013	179	13817	34679	18	2	822	810	988	983	70	154	11804	34500	-	-
560	116	TRUE	FALSE	1162	124	10094	17756	18	3	1196	1207	987	990	95	155	8932	17632	-	-
562	1073	TRUE	FALSE	34599	1880	94687	166049	98	12	542	558	1010	1002	154	165	60088	164169	-	-
563	1426	TRUE	TRUE	33521	4855	124785	223033	76	11	1104	1102	1020	999	140	164	91264	218178	-	-
564	1374	TRUE	FALSE	45084	2573	145386	211421	104	6	1197	1210	1018	1006	177	158	100302	208848	-	-
565	1767	TRUE	FALSE	68550	6998	158667	279116	83	6	471	445	996	994	134	160	90117	272118	-	-
566	469	TRUE	FALSE	9793	761	36995	71111	40	3	702	702	999	990	98	153	27202	70350	-	-
567	901	TRUE	FALSE	23562	1886	89335	138838	66	4	1182	1173	1024	1033	139	156	65773	136952	-	-
568	361	TRUE	FALSE	1746	666	20518	55177	11	3	854	851	1002	992	63	154	18772	54511	-	-
570	265	TRUE	FALSE	2139	430	16184	40975	14	4	571	575	1005	1005	67	157	14045	40545	-	-
571	519	TRUE	FALSE	11706	703	42327	78553	50	2	673	675	1005	999	109	152	30621	77850	-	-
572	909	TRUE	FALSE	37841	1680	85109	138030	129	6	723	751	1007	1017	181	156	47268	136350	-	-
573	1491	TRUE	FALSE	30030	1933	104580	231547	48	4	1064	1056	1013	1006	98	158	74550	229614	-	-
574	133	TRUE	FALSE	329	173	5649	21586	4	3	366	356	1006	1007	44	164	5320	21413	-	-
575	286	TRUE	FALSE	1940	317	15096	44361	12	5	1016	1020	1013	1012	58	159	13156	44044	-	-
576	888	TRUE	FALSE	18924	1511	59772	137375	68	3	974	975	1029	1006	114	156	40848	135864	-	-
577	221	TRUE	FALSE	1136	340	12186	33711	11	3	892	903	1018	1011	61	154	11050	33371	-	-
578	337	TRUE	FALSE	2452	260	17954	52158	13	2	1001	1002	1020	1008	59	156	15502	51898	-	-
580	101	TRUE	FALSE	374	70	5929	15321	6	2	871	871	1017	1013	61	153	5555	15251	-	-
581	742	TRUE	FALSE	8417	757	44775	112799	28	3	843	856	1034	1037	77	154	36358	112042	-	-
582	799	TRUE	FALSE	11806	878	54952	120727	41	6	669	660	1029	1023	95	156	43146	119849	-	-
583	892	TRUE	FALSE	26622	3734	68546	134858	95	6	720	719	1029	1018	142	153	41924	131124	-	-
584	1002	TRUE	TRUE	19677	2198	71781	155504	57	8	570	573	1036	1036	109	161	52104	153306	-	-
585	346	TRUE	FALSE	2106	450	21482	53042	12	6	624	614	1030	1023	68	158	19376	52592	-	-
586	261	TRUE	FALSE	2052	276	15363	39687	14	2	869	866	1028	1020	65	153	13311	39411	-	-
587	912	TRUE	FALSE	18726	821	58854	138517	61	4	890	898	1052	1055	105	155	40128	137696	-	-
588	386	TRUE	FALSE	4124	551	23810	58837	21	2	780	784	1022	1021	72	153	19686	58286	-	-
590	507	TRUE	FALSE	8571	1850	41019	79928	26	10	478	488	1031	1030	90	164	32448	78078	-	-
591	1005	TRUE	TRUE	26626	3081	88936	157848	60	9	511	509	1032	1032	122	163	62310	154767	-	-
594	2655	TRUE	FALSE	93257	3867	212732	404770	152	22	949	971	1087	1087	197	173	119475	400903	-	-
595	642	TRUE	FALSE	17355	1873	48171	98815	76	11	982	997	1050	1054	124	162	30816	96942	-	-
596	751	TRUE	FALSE	11346	1359	53402	116262	35	3	1210	1219	1053	1040	91	156	42056	114903	-	-
597	265	TRUE	FALSE	3406	206	15596	40221	31	1	936	933	1052	1037	77	152	12190	40015	-	-
598	438	TRUE	FALSE	5634	1289	27972	68303	28	6	1036	1029	1058	1050	79	159	22338	67014	-	-
599	1424	TRUE	FALSE	29867	3134	118155	222428	45	7	481	464	1055	1041	107	161	88288	219294	-	-
600	638	TRUE	FALSE	4192	1460	37368	98436	11	7	624	622	1050	1070	63	159	33176	96976	-	-
601	793	TRUE	FALSE	13651	1368	50129	121111	50	4	748	740	1054	1043	96	155	36478	119743	-	-
602	633	TRUE	FALSE	6924	1119	40473	97967	23	4	1073	1069	1049	1066	76	157	33549	96848	-	-
603	1015	TRUE	FALSE	14806	1562	40181	155840	42	4	1335	1329	1057	1061	67	156	25375	154278	-	-
604	402	TRUE	FALSE	2841	622	21735	61324	13	20	831	822	1053	1048	60	171	18894	60702	-	-
605	445	TRUE	FALSE	3260	858	25065	67608	14	3	693	687	1053	1047	63	153	21805	66750	-	-
606	215	TRUE	FALSE	2026	148	12991	32613	15	3	674	680	1051	1050	66	154	10965	32465	-	-
607	50	TRUE	FALSE	239	37	3189	7737	10	2	1068	1065	1049	1046	69	156	2950	7700	-	-
608	1053	TRUE	FALSE	15657	2489	77784	163594	29	6	1149	1165	1056	1070	88	159	62127	161105	-	-

609	520	TRUE	FALSE	8003	810	39723	80890	27	3	1198	1189	1055	1049	88	157	31720	80080	-	-
610	989	TRUE	FALSE	21985	1979	54622	152307	43	8	1267	1256	1061	1060	76	160	32637	150328	-	-
611	670	TRUE	FALSE	12359	1611	43849	102781	40	11	1020	1008	1067	1058	87	162	31490	101170	-	-
612	718	TRUE	FALSE	9675	691	44857	110545	20	3	557	552	1063	1064	69	156	35182	109854	-	-
613	927	TRUE	FALSE	40945	1483	82660	141460	158	5	980	996	1067	1056	203	156	41715	139977	-	-
614	664	TRUE	FALSE	7653	778	40189	101706	25	3	800	807	1072	1076	74	155	32536	100928	-	-
615	1018	TRUE	FALSE	34706	2207	81534	155925	128	6	890	907	1069	1080	174	157	46828	153718	-	-
616	297	TRUE	FALSE	5140	685	19693	45829	35	4	1038	1045	1068	1063	84	156	14553	45144	-	-
617	1394	TRUE	FALSE	26816	1738	88152	213626	57	4	856	859	1091	1080	101	156	61336	211888	-	-
618	252	TRUE	FALSE	1491	535	16107	39343	11	5	1107	1108	1071	1068	69	159	14616	38808	-	-
619	1280	TRUE	FALSE	16793	1545	106392	197385	48	5	488	512	1079	1092	118	158	89599	195840	-	-
620	505	TRUE	FALSE	6914	505	37719	78273	27	11	472	454	1085	1094	88	165	30805	77768	-	-
621	783	TRUE	FALSE	12257	1761	38879	122343	33	8	1264	1266	1074	1089	67	162	26622	120582	-	-
622	167	TRUE	FALSE	1621	280	11474	25998	16	3	1150	1151	1075	1082	75	157	9853	25718	-	-
623	195	TRUE	FALSE	1681	364	12796	30394	15	4	1170	1171	1082	1083	72	158	11115	30030	-	-
624	724	TRUE	FALSE	11161	1273	43017	110597	34	5	1032	1041	1088	1097	78	156	31856	109324	-	-
626	221	TRUE	FALSE	1786	369	10626	33961	13	3	890	892	1088	1088	53	155	8840	33592	-	-
627	297	TRUE	FALSE	2105	447	18440	45888	13	5	1108	1099	1095	1097	68	158	16335	45441	-	-
628	1004	TRUE	FALSE	20150	1759	70349	153363	55	6	635	650	1109	1112	105	157	50199	151604	-	-
629	1585	TRUE	FALSE	40687	4022	100917	251282	73	5	385	371	1114	1093	111	161	60230	247260	-	-
631	365	TRUE	FALSE	3237	321	24407	56531	17	6	1158	1150	1101	1105	75	160	21170	56210	-	-
632	593	TRUE	FALSE	5606	776	37628	90912	20	2	593	589	1110	1093	74	154	32022	90136	-	-
635	11	TRUE	FALSE	14	4	894	1698	2	1	472	474	1095	1095	82	155	880	1694	-	-
636	1857	TRUE	FALSE	79610	7510	178031	291631	131	10	471	469	1116	1128	184	163	98421	284121	-	-
637	940	TRUE	FALSE	25922	1623	70102	143562	78	7	1065	1055	1117	1118	125	158	44180	141939	-	-
638	1104	TRUE	FALSE	24674	1794	74354	167394	72	3	1031	1029	1122	1100	117	153	49680	165600	-	-
639	1346	TRUE	FALSE	43767	3435	129911	208027	115	7	497	483	1126	1129	179	159	86144	204592	-	-
640	588	TRUE	FALSE	9380	676	34664	89464	33	3	973	964	1116	1108	76	154	25284	88788	-	-
641	346	TRUE	FALSE	4215	549	23245	53141	25	7	613	606	1116	1110	80	159	19030	52592	-	-
642	465	TRUE	FALSE	5716	877	31756	72022	24	84	1160	1160	1124	1117	80	237	26040	71145	-	-
644	266	TRUE	FALSE	2181	266	15481	40432	17	3	961	948	1119	1114	67	154	13300	40166	-	-
645	340	TRUE	FALSE	2544	330	16824	51670	14	2	939	936	1120	1111	56	153	14280	51340	-	-
646	816	TRUE	FALSE	14659	1225	54643	126073	36	5	1213	1218	1129	1116	85	158	39984	124848	-	-
647	858	TRUE	FALSE	14748	1088	64512	131504	39	5	535	526	1128	1119	97	157	49764	130416	-	-
648	1099	TRUE	FALSE	29840	1410	77097	168458	91	12	712	707	1128	1118	134	164	47257	167048	-	-
649	389	TRUE	FALSE	3209	698	24604	59826	14	4	1113	1112	1126	1118	69	156	21395	59128	-	-
652	284	TRUE	FALSE	3268	528	17752	44548	21	3	426	417	1131	1123	72	158	14484	44020	-	-
653	662	TRUE	FALSE	11955	2113	43731	101373	37	8	626	624	1127	1131	85	158	31776	99260	-	-
655	583	TRUE	FALSE	11226	467	33963	88500	40	2	1008	1005	1136	1128	79	153	22737	88033	-	-
656	619	TRUE	FALSE	7069	848	41114	94936	22	2	1140	1135	1136	1117	77	154	34045	94088	-	-
657	161	TRUE	FALSE	1079	182	11061	24815	11	2	1197	1192	1128	1118	73	155	9982	24633	-	-
658	472	TRUE	FALSE	8818	386	26282	71655	35	3	995	988	1135	1142	72	154	17464	71269	-	-
659	671	TRUE	FALSE	20002	695	56907	102016	87	3	1076	1090	1134	1125	142	154	36905	101321	-	-
660	506	TRUE	FALSE	3713	555	23953	76961	13	3	831	818	1131	1133	53	154	20240	76406	-	-
661	306	TRUE	FALSE	1242	513	14094	46719	7	3	902	903	1128	1132	49	154	12852	46206	-	-
662	665	TRUE	FALSE	13283	1160	49193	102235	58	3	1164	1161	1138	1129	112	155	35910	101075	-	-
663	291	TRUE	FALSE	4143	233	17820	44174	27	2	971	980	1129	1125	74	153	13677	43941	-	-
664	315	TRUE	FALSE	1885	608	16690	48173	11	9	659	660	1131	1129	58	160	14805	47565	-	-
665	503	TRUE	FALSE	3999	847	24119	76800	14	3	800	813	1136	1132	54	154	20120	75953	-	-
666	222	TRUE	FALSE	3891	219	14103	33519	35	2	1053	1049	1136	1129	81	152	10212	33300	-	-
667	103	TRUE	FALSE	1416	169	7596	15619	25	2	1059	1051	1135	1128	85	152	6180	15450	-	-
668	399	TRUE	FALSE	4139	378	27680	62622	17	3	426	431	1143	1133	76	159	23541	62244	-	-
669	216	TRUE	FALSE	1958	171	13838	33003	16	3	1118	1122	1138	1136	71	155	11880	32832	-	-
670	333	TRUE	FALSE	4884	1048	24864	51997	29	7	472	473	1140	1137	89	160	19980	50949	-	-
671	1288	TRUE	FALSE	32601	1308	90560	197083	79	7	585	593	1156	1173	124	159	57959	195775	-	-
672	804	TRUE	FALSE	11913	801	41661	122204	38	4	825	824	1158	1150	75	155	29748	121403	-	-
673	1520	TRUE	FALSE	20437	3845	79717	244005	34	6	355	325	1162	1158	73	164	59280	240160	-	-
674	1312	TRUE	FALSE	30719	1267	93695	200684	64	10	540	553	1162	1156	112	162	62976	199417	-	-
675	420	TRUE	FALSE	8126	581	27866	63581	39	3	1058	1059	1156	1150	86	153	19740	63000	-	-
676	1176	TRUE	FALSE	29820	2698	66276	192034	84	5	285	272	1170	1180	115	166	36456	189336	-	-

677	384	TRUE	FALSE	5885	684	27005	58284	41	4	1090	1095	1159	1154	96	154	21120	57600	-	-
678	319	TRUE	FALSE	1936	252	14377	48421	11	2	977	981	1155	1149	50	153	12441	48169	-	-
679	757	TRUE	FALSE	10495	1737	46831	119072	26	4	424	409	1160	1154	74	159	36336	117335	-	-
680	306	TRUE	FALSE	3342	456	16806	46662	19	5	734	740	1160	1161	63	156	13464	46206	-	-
681	204	TRUE	FALSE	2447	272	14891	30872	20	3	1075	1080	1158	1151	81	153	12444	30600	-	-
682	1329	TRUE	FALSE	31788	1710	82290	203709	95	4	712	716	1174	1157	133	156	50502	201999	-	-
683	718	TRUE	FALSE	9931	587	37215	109723	33	6	649	641	1172	1176	71	158	27284	109136	-	-
684	187	TRUE	FALSE	2997	218	13843	28642	35	2	1146	1149	1172	1154	93	154	10846	28424	-	-
685	1427	TRUE	FALSE	43013	2072	120071	218976	111	12	1167	1178	1181	1179	165	164	77058	216904	-	-
686	1414	TRUE	FALSE	7703	2823	35983	216337	10	10	1345	1333	1177	1187	30	161	28280	213514	-	-
687	944	TRUE	FALSE	23065	1012	67433	143552	58	3	1114	1114	1181	1164	105	154	44368	142540	-	-
688	1678	TRUE	FALSE	74937	2768	152125	256146	169	11	1132	1143	1187	1197	215	162	77188	253378	-	-
689	33	TRUE	FALSE	215	52	1964	5068	12	2	603	605	1162	1162	65	154	1749	5016	-	-
690	714	TRUE	FALSE	5422	1856	45406	111098	13	14	474	485	1167	1197	69	167	39984	109242	-	-
691	526	TRUE	FALSE	6046	797	24456	79697	22	3	873	882	1173	1170	57	153	18410	78900	-	-
692	338	TRUE	FALSE	3699	438	14853	54856	17	4	317	318	1173	1183	50	165	11154	54418	-	-
693	592	TRUE	FALSE	7848	798	32712	90782	32	3	617	625	1175	1168	74	155	24864	89984	-	-
694	1147	TRUE	FALSE	19257	1100	80048	176591	47	6	1205	1201	1184	1169	100	159	60791	175491	-	-
695	1577	TRUE	FALSE	52899	1514	112825	238058	99	5	1061	1090	1187	1180	137	155	59926	236544	-	-
696	464	TRUE	FALSE	5110	459	19958	70059	20	12	981	980	1186	1178	52	162	14848	69600	-	-
697	534	TRUE	FALSE	5952	1078	29448	82246	26	4	589	579	1188	1197	70	156	23496	81168	-	-
698	289	TRUE	FALSE	1803	590	16253	44807	11	9	537	548	1193	1198	61	162	14450	44217	-	-
699	1753	TRUE	FALSE	53825	5936	104662	270639	100	9	916	917	1223	1243	129	160	50837	264703	-	-
701	1169	TRUE	FALSE	7673	3761	67292	187294	11	7	468	476	1217	1235	62	164	59619	183533	-	-
702	1108	TRUE	FALSE	12992	743	47340	170238	20	2	1240	1258	1210	1205	51	155	34348	169495	-	-
704	939	TRUE	FALSE	6855	2025	40659	145692	11	7	717	710	1216	1238	47	160	33804	143667	-	-
705	408	TRUE	FALSE	2684	342	16556	61542	13	3	1026	1022	1220	1214	47	153	13872	61200	-	-
706	830	TRUE	FALSE	10741	1476	46431	127636	30	7	1197	1208	1225	1212	73	159	35690	126160	-	-
707	532	TRUE	FALSE	4895	1307	25643	88555	18	6	380	390	1229	1229	57	170	20748	87248	-	-
708	174	TRUE	TRUE	1384	739	9214	27013	12	8	1136	1130	1224	1224	57	159	7830	26274	-	-
709	530	TRUE	FALSE	6296	1429	21136	89409	29	9	286	288	1235	1226	57	175	14840	87980	-	-
710	636	TRUE	FALSE	5682	995	32394	99998	18	5	594	596	1242	1255	60	161	26712	99003	-	-
711	556	TRUE	FALSE	4704	1394	28056	88128	16	5	553	542	1238	1239	58	161	23352	86734	-	-
712	1729	TRUE	FALSE	22307	3891	91467	283989	29	32	409	394	1250	1255	69	194	69160	280098	-	-
713	305	TRUE	FALSE	2497	248	15612	46303	17	3	1137	1130	1241	1232	60	154	13115	46055	-	-
714	341	TRUE	FALSE	1877	457	16540	51948	11	5	1160	1163	1243	1247	54	156	14663	51491	-	-
715	137	TRUE	FALSE	504	396	4751	23001	5	9	316	310	1236	1241	36	174	4247	22605	-	-
716	325	TRUE	FALSE	2663	569	16313	49969	14	3	1200	1203	1255	1252	56	155	13650	49400	-	-
717	430	TRUE	FALSE	2154	865	16774	64935	10	6	1062	1054	1260	1256	44	155	14620	64070	-	-
718	678	TRUE	FALSE	8120	1221	35240	102921	32	13	1138	1135	1265	1280	72	163	27120	101700	-	-
719	349	TRUE	FALSE	1744	616	13959	52617	10	15	1090	1089	1263	1255	45	164	12215	52001	-	-
720	147	TRUE	FALSE	420	270	5712	23202	5	4	857	855	1263	1266	41	160	5292	22932	-	-
721	1083	TRUE	TRUE	20270	23766	61424	197046	43	58	617	620	1283	1282	81	218	41154	173280	-	-
722	846	TRUE	TRUE	13882	7107	50260	142467	37	28	594	607	1294	1284	80	188	36378	135360	-	-
723	288	TRUE	FALSE	1846	521	13654	44297	13	7	1197	1205	1271	1270	54	159	11808	43776	-	-
724	1405	TRUE	TRUE	29227	41270	78402	263260	61	97	709	709	1283	1281	96	255	49175	221990	-	-
725	512	TRUE	FALSE	4521	1756	22441	82140	19	7	858	851	1283	1277	54	164	17920	80384	-	-
726	914	TRUE	TRUE	7224	7457	50182	154611	14	21	553	557	1288	1293	61	182	42958	147154	-	-
727	735	TRUE	FALSE	3984	1866	27504	117261	14	10	804	797	1284	1292	46	167	23520	115395	-	-
728	201	TRUE	FALSE	1109	174	5129	31329	9	2	1340	1348	1280	1275	29	157	4020	31155	-	-
729	658	TRUE	FALSE	3297	545	28301	101218	10	5	1216	1223	1293	1278	48	158	25004	100673	-	-
730	1065	TRUE	FALSE	8927	2116	32357	167191	15	5	1343	1343	1303	1303	37	160	23430	165075	-	-
731	173	TRUE	FALSE	878	317	6760	27824	10	4	858	852	1297	1301	44	163	5882	27507	-	-
732	1230	TRUE	TRUE	30503	20332	82163	220822	55	45	598	614	1307	1304	97	208	51660	200490	-	-
733	1139	TRUE	TRUE	19287	27805	60291	211184	42	63	701	696	1310	1309	78	224	41004	183379	-	-
734	426	TRUE	TRUE	2892	2293	25044	73435	12	12	542	546	1309	1319	64	179	22152	71142	-	-
735	57	TRUE	FALSE	152	61	2204	9238	5	2	857	856	1300	1303	41	163	2052	9177	-	-
736	543	TRUE	FALSE	3648	794	18309	84416	11	3	1267	1280	1320	1309	38	157	14661	83622	-	-
737	1553	TRUE	FALSE	29677	2787	93350	235737	45	4	1140	1147	1350	1322	86	154	63673	232950	-	-
738	220	TRUE	FALSE	1294	365	8334	35125	13	3	887	882	1328	1320	45	161	7040	34760	-	-

740	711	TRUE	TRUE	7721	2394	33317	114732	24	14	681	682	1330	1319	60	172	25596	112338	-	-
741	1601	TRUE	TRUE	47687	25713	113328	283474	88	46	594	585	1336	1351	129	207	65641	257761	-	-
742	960	TRUE	TRUE	10170	8577	56250	163137	22	21	549	538	1340	1354	70	182	46080	154560	-	-
743	2280	TRUE	FALSE	60554	6843	142634	348843	75	17	1215	1212	1354	1345	111	167	82080	342000	-	-
745	666	TRUE	FALSE	9362	787	30008	106013	32	7	888	905	1340	1344	63	165	20646	105226	-	-
747	703	TRUE	FALSE	13713	12006	34803	136437	38	78	284	275	1354	1349	68	255	21090	124431	-	-
748	1174	TRUE	FALSE	18957	8589	55351	204647	44	15	788	779	1357	1345	75	182	36394	196058	-	-
749	617	TRUE	FALSE	8593	7989	27103	112262	28	86	297	297	1353	1357	58	255	18510	104273	-	-
750	238	TRUE	TRUE	1576	2328	13000	41598	14	17	535	536	1351	1355	62	182	11424	39270	-	-
751	723	TRUE	FALSE	9332	1203	32468	115434	29	4	885	872	1360	1364	61	162	23136	114231	-	-
752	207	TRUE	FALSE	914	138	8159	33042	7	2	695	702	1353	1351	42	161	7245	32904	-	-
753	526	TRUE	FALSE	4288	594	21120	84228	15	5	723	726	1362	1355	47	164	16832	83634	-	-
754	527	TRUE	FALSE	3674	2683	22119	90163	15	33	375	372	1361	1353	50	199	18445	87480	-	-
755	1250	TRUE	FALSE	24014	2287	67764	193535	50	7	1070	1085	1376	1383	85	160	43750	191248	-	-
756	1464	TRUE	FALSE	50919	5387	109479	223523	127	9	1134	1125	1379	1372	167	158	58560	218136	-	-
757	223	TRUE	FALSE	912	278	7825	34173	6	3	1264	1265	1360	1367	37	155	6913	33895	-	-
758	396	TRUE	FALSE	1387	405	13267	61386	6	2	994	992	1370	1363	36	156	11880	60981	-	-
760	818	TRUE	FALSE	5655	4738	36739	140524	14	19	379	375	1386	1397	52	185	31084	135786	-	-
762	216	TRUE	FALSE	1036	291	8164	34419	8	2	879	880	1378	1371	41	160	7128	34128	-	-
763	393	TRUE	FALSE	3419	658	19925	59215	18	3	1190	1198	1381	1381	60	152	16506	58557	-	-
764	555	TRUE	FALSE	6964	3132	29164	91932	25	10	414	410	1389	1391	65	170	22200	88800	-	-
765	223	TRUE	FALSE	1042	287	8401	34406	8	2	1031	1026	1391	1385	41	155	7359	34119	-	-
766	333	TRUE	FALSE	796	467	8122	50750	4	3	1345	1342	1407	1391	26	154	7326	50283	-	-
767	804	TRUE	FALSE	6288	4331	42468	130559	16	22	489	496	1407	1423	61	179	36180	126228	-	-
768	1046	TRUE	FALSE	13126	1615	60196	158515	27	6	1196	1188	1442	1436	72	156	47070	156900	-	-
769	851	TRUE	FALSE	7735	730	45179	130076	16	12	1159	1163	1448	1417	60	164	37444	129346	-	-
770	428	TRUE	FALSE	2857	349	22117	65832	12	2	1122	1123	1432	1429	57	155	19260	65483	-	-
771	414	TRUE	FALSE	2783	308	19343	63650	10	2	1073	1080	1434	1423	50	155	16560	63342	-	-
772	252	TRUE	FALSE	1670	392	11246	38696	10	3	1047	1040	1435	1424	48	155	9576	38304	-	-
773	209	TRUE	FALSE	308	441	5533	31791	3	12	1343	1342	1436	1434	28	162	5225	31350	-	-
774	1274	TRUE	FALSE	24456	2133	79238	200852	54	16	526	511	1445	1434	97	172	54782	198719	-	-
775	196	TRUE	FALSE	244	311	7300	30495	2	4	862	876	1431	1429	38	158	7056	30184	-	-
777	1462	TRUE	FALSE	14541	4145	74483	233679	19	11	460	465	1446	1431	60	168	59942	229534	-	-
778	804	TRUE	FALSE	7261	1157	29773	121757	16	16	1262	1281	1447	1448	44	166	22512	120600	-	-
779	577	TRUE	FALSE	8141	831	31221	88511	31	4	1045	1046	1459	1449	71	156	23080	87680	-	-
780	68	TRUE	FALSE	81	129	2597	10533	2	4	862	865	1451	1454	39	157	2516	10404	-	-
781	279	TRUE	FALSE	3460	218	18247	42347	19	2	1195	1193	1451	1453	72	153	14787	42129	-	-
782	1448	TRUE	FALSE	25833	6455	88097	235238	48	97	409	420	1468	1472	91	255	62264	228783	-	-
784	128	TRUE	FALSE	320	186	4928	19770	4	2	792	789	1460	1458	40	155	4608	19584	-	-
785	253	TRUE	FALSE	1442	250	9285	39212	9	4	658	652	1467	1464	40	158	7843	38962	-	-
797	1991	TRUE	FALSE	157102	8245	252668	296940	207	11	576	573	433	456	255	156	95566	288695	-	-
799	545	TRUE	FALSE	60760	2656	88555	82226	204	10	562	555	436	452	255	156	27795	79570	-	-
800	1133	TRUE	FALSE	40205	6729	103653	173280	108	11	543	506	437	435	164	158	63448	166551	-	-
801	465	TRUE	FALSE	33536	2278	57716	70168	160	11	554	552	437	462	212	157	24180	67890	-	-
802	261	TRUE	FALSE	4715	499	22463	38344	34	4	635	626	436	435	102	149	17748	37845	-	-
803	465	TRUE	FALSE	10051	794	35626	69614	46	4	634	628	429	415	101	152	25575	68820	-	-
804	1880	TRUE	FALSE	221460	12006	336140	282726	194	13	786	780	443	424	255	157	114680	270720	-	-
805	2097	TRUE	FALSE	233864	12085	355490	314053	197	11	815	828	441	443	255	155	121626	301968	-	-
806	160	TRUE	FALSE	536	530	5816	25330	7	6	1045	1045	344	341	40	161	5280	24800	-	-
807	136	TRUE	FALSE	526	494	4606	21846	9	17	1072	1077	346	348	39	174	4080	21352	-	-
808	114	TRUE	FALSE	358	275	4006	17603	6	12	1152	1152	344	346	38	164	3648	17328	-	-
809	104	TRUE	FALSE	322	199	3442	15903	7	5	1174	1179	345	344	37	156	3120	15704	-	-
810	120	TRUE	FALSE	492	192	4452	18912	8	4	1018	1016	344	341	41	160	3960	18720	-	-
811	500	TRUE	FALSE	3708	1259	21707	77259	39	6	1021	1014	317	324	75	158	17999	76000	-	-
812	480	TRUE	FALSE	3527	1421	21767	74861	24	8	990	981	317	322	62	161	18240	73440	-	-
817	124	TRUE	FALSE	337	817	4181	19913	5	24	1075	1076	357	354	36	178	3844	19096	-	-
818	124	TRUE	FALSE	250	937	4342	19909	4	36	1096	1099	359	355	37	189	4092	18972	-	-
819	142	TRUE	FALSE	419	550	5105	22418	5	10	1051	1047	355	353	38	164	4686	21868	-	-
820	124	TRUE	FALSE	594	353	4686	19449	8	11	1019	1021	357	353	41	165	4092	19096	-	-
821	124	TRUE	FALSE	367	485	4707	19333	9	7	1155	1156	324	326	44	159	4340	18848	-	-

823	90	TRUE	FALSE	539	88	3689	13768	12	4	1155	1156	334	329	47	156	3150	13680	-	-
824	157	TRUE	FALSE	1266	385	6447	24244	17	15	1162	1168	331	331	50	167	5181	23859	-	-
842	98	TRUE	FALSE	464	322	3306	15316	10	10	1094	1089	297	298	39	163	2842	14994	-	-
843	63	TRUE	FALSE	376	195	2203	9834	10	6	1109	1109	298	298	39	159	1827	9639	-	-
846	104	TRUE	FALSE	820	574	4252	16486	20	10	1095	1089	304	303	53	163	3432	15912	-	-
847	112	TRUE	FALSE	1181	682	4877	17706	29	10	1112	1110	305	304	62	162	3696	17024	-	-
848	79	TRUE	FALSE	519	121	3126	12129	12	4	1094	1092	310	312	45	156	2607	12008	-	-
850	93	TRUE	FALSE	727	100	3889	14236	14	2	1110	1106	311	309	48	154	3162	14136	-	-
852	100	TRUE	FALSE	844	82	4344	15282	17	1	1127	1121	311	309	52	153	3500	15200	-	-
853	190	TRUE	FALSE	2044	401	8124	29091	23	6	1147	1153	312	303	55	157	6080	28690	-	-
854	178	TRUE	FALSE	2131	540	7471	27596	38	7	1127	1122	306	304	68	159	5340	27056	-	-
855	140	TRUE	FALSE	1434	441	5774	21721	30	6	1142	1147	306	303	61	158	4340	21280	-	-
856	120	TRUE	FALSE	545	715	4385	19075	11	12	1076	1078	304	300	43	165	3840	18360	-	-
857	64	TRUE	FALSE	150	99	2454	9827	5	2	1072	1067	309	308	41	154	2304	9728	-	-
858	124	TRUE	FALSE	559	250	4279	19222	7	6	1050	1047	292	297	37	159	3720	18972	-	-
859	70	TRUE	FALSE	193	130	2293	10770	5	4	1031	1034	291	295	35	156	2100	10640	-	-
860	171	TRUE	FALSE	809	220	4913	26383	8	3	1076	1081	285	283	32	156	4104	26163	-	-
861	153	TRUE	FALSE	655	119	4480	23220	8	5	1153	1157	289	285	33	156	3825	23101	-	-
862	148	TRUE	FALSE	579	440	5759	23380	7	71	1002	1005	331	336	42	226	5180	22940	-	-
863	151	TRUE	FALSE	647	150	5630	23404	8	3	1031	1029	329	337	41	157	4983	23254	-	-
864	181	TRUE	FALSE	608	281	7305	28336	6	4	974	971	330	329	43	159	6697	28055	-	-
870	1208	TRUE	FALSE	10702	3532	63854	193185	18	24	863	908	315	300	62	181	53152	189653	-	-
871	75	TRUE	FALSE	398	76	3698	11851	10	2	867	872	306	303	54	159	3300	11775	-	-
874	111	TRUE	FALSE	673	181	6001	17497	11	4	884	875	333	336	59	160	5328	17316	-	-
876	110	TRUE	FALSE	580	190	5970	17350	10	4	907	906	334	331	59	160	5390	17160	-	-
877	135	TRUE	FALSE	1246	136	7591	21331	17	3	838	846	316	314	64	160	6345	21195	-	-
879	5103	TRUE	FALSE	89002	23420	308431	799076	34	44	824	723	329	337	77	196	219429	775656	-	-
880	131	TRUE	FALSE	399	250	4198	20162	7	4	1191	1194	374	379	36	156	3799	19912	-	-
881	184	TRUE	FALSE	466	713	5986	28681	5	19	1174	1183	365	360	35	171	5520	27968	-	-
882	142	TRUE	FALSE	715	108	5827	21265	9	1	1142	1140	436	431	45	150	5112	21157	-	-
883	263	TRUE	FALSE	845	777	10576	39175	7	5	1148	1141	447	445	44	151	9731	38398	-	-
884	216	TRUE	FALSE	701	613	6749	32581	6	8	1199	1200	411	405	34	156	6048	31968	-	-
885	368	TRUE	TRUE	8431	2858	19839	57690	46	44	1430	1432	369	360	77	193	11408	54832	-	-
887	1000	TRUE	TRUE	22446	4259	71446	153259	89	24	939	940	436	439	138	173	49000	149000	-	-
888	661	TRUE	TRUE	23401	3342	57112	101170	95	15	905	909	437	437	146	163	33711	97828	-	-
889	392	TRUE	FALSE	8354	610	29914	58626	69	4	903	920	446	442	124	152	21560	58016	-	-
890	272	TRUE	FALSE	1833	343	11625	41143	11	8	1014	1018	423	420	47	158	9792	40800	-	-
891	183	TRUE	FALSE	878	293	7283	27377	8	4	1076	1069	431	433	43	152	6405	27084	-	-
892	720	TRUE	FALSE	23926	939	72166	103177	65	13	810	793	506	509	132	155	48240	102238	-	-
893	1186	TRUE	FALSE	45896	2500	114684	169726	114	5	818	827	506	482	172	146	68788	167226	-	-
894	1125	TRUE	FALSE	32899	3057	89149	163932	89	6	931	928	501	479	139	149	56250	160875	-	-
895	1737	TRUE	FALSE	81942	4182	163581	256047	132	6	944	969	502	505	179	151	81639	251865	-	-
896	370	TRUE	FALSE	3412	1831	21912	56591	22	12	996	980	517	511	72	160	99500	54760	-	-
897	2021	TRUE	FALSE	89241	7300	178165	302306	137	18	1011	1021	505	510	181	164	88924	295006	-	-
898	312	TRUE	FALSE	2510	1240	18110	47104	16	21	912	919	457	460	66	168	15600	45864	-	-
899	216	TRUE	FALSE	1006	971	11158	32291	9	12	594	589	478	479	56	157	10152	31320	-	-
900	585	TRUE	FALSE	2709	443	34299	88192	10	2	631	642	551	541	64	152	31590	87749	-	-
901	570	TRUE	FALSE	1684	957	31879	85887	6	5	608	616	544	545	59	154	30195	84930	-	-
902	294	TRUE	FALSE	1440	536	16434	44342	9	3	590	583	559	556	60	152	14994	43806	-	-
903	384	TRUE	FALSE	3026	540	24146	57756	14	3	658	650	554	557	69	152	21120	57216	-	-
904	253	TRUE	FALSE	1715	327	16136	38024	12	3	672	670	550	542	69	152	14421	37697	-	-
905	269	TRUE	FALSE	2320	934	20343	41822	14	8	546	552	602	604	81	160	18023	40888	-	-
907	86	TRUE	FALSE	598	134	6274	13120	12	2	561	558	583	577	78	153	5676	12986	-	-
908	444	TRUE	FALSE	8395	1915	37699	69403	37	13	563	565	598	599	103	165	29304	67488	-	-
910	775	TRUE	FALSE	40048	1713	84998	117188	137	5	680	668	586	594	195	154	44950	115475	-	-
911	1759	TRUE	FALSE	117999	3828	220021	267678	197	7	704	706	582	599	255	157	102022	263850	-	-
914	856	TRUE	FALSE	14183	1165	61263	129513	36	9	628	629	579	573	91	159	47080	128348	-	-
915	1105	TRUE	FALSE	38049	2842	98824	167487	83	6	666	658	586	599	138	155	60775	164645	-	-
916	313	TRUE	FALSE	1537	410	20317	47673	10	3	751	749	550	541	70	154	18780	47263	-	-
917	372	TRUE	FALSE	1872	680	24192	56852	9	4	805	806	591	595	69	155	22320	56172	-	-

918	250	TRUE	FALSE	998	255	12998	38245	9	3	974	983	556	557	57	155	12000	37990	-	-
919	307	TRUE	FALSE	2168	281	17518	47558	13	4	1030	1031	593	584	63	158	15350	47277	-	-
920	572	TRUE	FALSE	6122	928	20422	86156	28	4	1330	1331	566	578	53	153	14300	85228	-	-
921	930	TRUE	FALSE	10841	2812	30371	144172	23	20	1335	1336	601	612	44	172	19530	141360	-	-
922	635	TRUE	FALSE	7271	1229	28226	97114	28	6	1259	1249	544	551	61	157	20955	95885	-	-
924	864	TRUE	FALSE	10884	2847	37668	133311	22	22	1244	1267	560	571	53	173	26784	130464	-	-
925	368	TRUE	FALSE	1206	387	19606	57059	7	3	834	831	691	694	57	157	18400	56672	-	-
926	417	TRUE	FALSE	2149	420	22165	64638	10	11	858	864	693	688	58	165	20016	64218	-	-
927	149	TRUE	FALSE	815	1898	8414	25142	9	23	834	832	668	670	60	179	7599	23244	-	-
928	147	TRUE	FALSE	950	1591	8300	24376	10	21	840	847	664	671	60	176	7350	22785	-	-
929	366	TRUE	TRUE	2091	3630	25149	59628	12	24	629	626	672	675	75	177	23058	55998	-	-
930	620	TRUE	FALSE	4313	446	35933	94685	13	3	686	685	813	805	64	155	31620	94239	-	-
931	181	TRUE	FALSE	759	178	9990	27690	11	2	673	663	811	809	62	154	9231	27512	-	-
932	348	TRUE	FALSE	3531	436	20235	52636	18	3	843	844	786	788	66	153	16704	52200	-	-
933	331	TRUE	FALSE	2465	289	19346	50270	17	2	849	848	779	770	68	153	16881	49981	-	-
934	306	TRUE	FALSE	1332	490	17244	46696	8	4	799	795	772	766	60	155	15912	46206	-	-
935	587	TRUE	FALSE	6023	1227	17763	90451	18	4	1332	1333	840	835	38	156	11740	89224	-	-
936	456	TRUE	FALSE	4275	954	13851	70266	16	3	1338	1331	857	850	37	155	9576	69312	-	-
937	555	TRUE	FALSE	17642	1538	60932	87008	85	5	1189	1191	844	827	163	159	43290	85470	-	-
939	285	TRUE	FALSE	13837	633	34357	44523	95	5	1189	1182	851	854	167	159	20520	43890	-	-
941	320	TRUE	FALSE	14963	337	35763	49937	92	3	1189	1172	859	862	157	158	20800	49600	-	-
942	692	TRUE	FALSE	18999	1175	59135	107743	81	4	1191	1191	870	880	139	158	40136	106568	-	-
943	462	TRUE	FALSE	7920	748	33330	70048	33	3	1098	1106	865	863	88	153	25410	69300	-	-
944	453	TRUE	TRUE	14067	1033	41247	68983	77	9	1074	1074	832	823	137	159	27180	67950	-	-
946	141	TRUE	FALSE	1757	148	8384	21298	21	2	1080	1069	856	853	68	152	6627	21150	-	-
947	104	TRUE	FALSE	1026	135	6850	15735	16	3	1083	1081	845	850	72	153	5824	15600	-	-
948	1539	TRUE	FALSE	112482	5358	169425	245442	203	6	1240	1227	814	803	240	162	56943	240084	-	-
949	1721	TRUE	FALSE	58001	7106	125120	275582	123	29	1239	1229	802	780	162	185	67119	268476	-	-
950	629	TRUE	FALSE	3755	1142	34576	93605	14	5	911	910	817	826	63	152	30821	92463	-	-
951	615	TRUE	FALSE	4285	1007	34420	93257	16	6	796	792	862	868	65	156	30135	92250	-	-
952	336	TRUE	FALSE	1321	259	18457	51330	7	7	650	655	823	825	58	159	17136	51071	-	-
953	681	TRUE	FALSE	5979	1410	29814	114456	20	11	288	282	964	955	55	177	23835	113046	-	-
954	372	TRUE	FALSE	8703	527	35487	58559	42	4	505	496	1013	1012	114	160	26784	58032	-	-
955	948	TRUE	FALSE	37417	1952	97141	148892	87	16	520	508	1012	989	150	171	59724	146940	-	-
956	129	TRUE	FALSE	773	198	8126	19935	12	2	610	616	969	968	69	155	7353	19737	-	-
957	264	TRUE	FALSE	1171	412	16747	40804	9	4	631	639	966	963	68	157	15576	40392	-	-
958	732	TRUE	FALSE	7871	682	42275	111214	21	2	880	889	986	960	68	153	34404	110532	-	-
959	476	TRUE	FALSE	9093	1301	41461	73653	38	6	1133	1132	1016	1011	106	158	32368	72352	-	-
961	370	TRUE	FALSE	7659	708	38369	56948	72	14	1145	1140	1011	1008	155	166	30710	56240	-	-
963	674	TRUE	TRUE	46238	1874	97462	104996	179	12	1145	1151	1018	1021	255	165	51224	103122	-	-
964	606	TRUE	FALSE	15379	1183	54769	93901	80	5	1149	1139	1032	1044	145	158	39390	92718	-	-
965	373	TRUE	TRUE	6719	3990	30964	60686	34	36	1215	1217	977	972	99	188	24245	56696	-	-
966	347	TRUE	TRUE	3560	1671	26115	54762	18	13	1214	1218	994	1003	83	166	22555	53091	-	-
968	2930	TRUE	FALSE	97978	6408	174158	448838	149	9	1268	1251	1035	1026	175	160	76180	442430	-	-
969	753	TRUE	FALSE	6131	1414	20438	115870	14	4	1347	1352	1127	1114	33	156	14307	114456	-	-
972	788	TRUE	FALSE	11614	1329	39982	121893	30	15	1252	1270	1130	1125	66	168	28368	120564	-	-
973	443	TRUE	FALSE	7672	895	26721	68674	41	3	1237	1233	1123	1107	84	156	19049	67779	-	-
974	410	TRUE	FALSE	9906	371	27126	63101	53	3	1238	1232	1131	1127	95	156	17220	62730	-	-
975	545	TRUE	FALSE	9113	1127	35818	83422	45	4	794	815	1094	1090	94	155	26705	82295	-	-
976	869	TRUE	FALSE	25811	1567	61440	132786	102	5	793	814	1104	1105	143	156	35629	131219	-	-
978	914	TRUE	FALSE	7961	800	54575	139728	21	3	695	719	1103	1099	72	155	46614	138928	-	-
979	1285	TRUE	FALSE	20443	3183	65418	198503	38	7	805	809	1221	1216	73	159	44975	195320	-	-
980	926	TRUE	FALSE	8279	1240	36985	144770	15	11	1042	1045	1350	1348	46	166	28706	143530	-	-
981	471	TRUE	FALSE	2084	2450	16214	79694	9	12	790	774	1388	1387	39	176	14130	77244	-	-
982	1547	TRUE	FALSE	34237	8632	94570	253058	44	19	441	446	1388	1392	83	177	60333	244426	-	-
983	526	TRUE	FALSE	5308	2216	26874	95318	17	31	257	243	1354	1354	58	208	21566	93102	-	-
984	596	TRUE	FALSE	9885	5642	22997	111134	28	55	242	226	1350	1344	50	232	13112	105492	-	-
985	3528	TRUE	FALSE	406444	9811	484060	535483	219	16	1199	1211	1521	1541	241	165	77616	525672	-	-
986	3071	TRUE	FALSE	377505	8581	405144	469231	177	9	1163	1160	1529	1500	186	159	27639	460650	-	-
987	3408	TRUE	FALSE	367801	12257	371209	523457	161	14	1120	1133	1536	1568	162	164	3408	511200	-	-

988	1906	TRUE	FALSE	146055	7074	163209	289162	96	17	1240	1253	1538	1516	105	165	17154	282088	-	-
990	1810	TRUE	FALSE	168021	6679	182501	279989	121	8	1079	1077	1528	1533	129	159	14480	273310	-	-
991	830	TRUE	FALSE	27209	2190	86969	127519	70	8	1086	1077	1524	1515	142	159	59760	125329	-	-
992	444	TRUE	FALSE	3300	671	22392	68603	16	5	843	850	1512	1507	59	158	19092	67932	-	-
993	558	TRUE	FALSE	4204	874	24292	86806	15	4	733	733	1510	1519	51	158	20088	85932	-	-
994	510	TRUE	FALSE	6259	378	26149	79934	22	2	611	608	1548	1535	61	158	19890	79556	-	-
996	468	TRUE	FALSE	4588	846	26584	73386	19	5	501	502	1496	1487	66	160	21996	72540	-	-
997	438	TRUE	FALSE	1922	661	18128	68989	8	3	709	711	1523	1525	45	159	16206	68328	-	-
998	786	TRUE	FALSE	6301	2032	40885	124648	17	6	455	477	1499	1496	61	162	34584	122616	-	-
1000	581	TRUE	FALSE	6809	937	34697	92735	24	5	413	420	1495	1504	72	163	27888	91798	-	-
1001	782	TRUE	FALSE	16333	1672	53087	126792	35	4	386	398	1540	1531	82	164	36754	125120	-	-
1002	590	TRUE	FALSE	4392	1051	23272	95451	13	37	333	323	1473	1475	45	197	18880	94400	-	-
1003	653	TRUE	FALSE	32018	956	70545	100212	114	4	722	726	994	983	173	156	38527	99256	-	-
1005	328	TRUE	TRUE	3382	1048	20438	50904	27	10	1016	1018	789	794	79	162	17056	49856	-	-
1006	834	TRUE	TRUE	15816	2400	55014	129168	49	8	1016	1020	783	781	96	160	39198	126768	-	-
1007	878	FALSE	TRUE	1272	4334	24978	140423	4	22	904	921	1242	1252	31	177	23706	136089	-	-
1008	1083	FALSE	TRUE	11846	5328	67079	168861	50	15	728	740	983	979	101	166	55233	163533	-	-
1009	229	FALSE	TRUE	150	858	13413	36124	4	21	1161	1169	899	890	62	175	13263	35266	-	-
1010	175	FALSE	TRUE	843	515	10293	27640	11	7	1105	1104	776	767	65	162	9450	27125	-	-
1011	382	TRUE	TRUE	5134	5123	26144	65861	36	49	1151	1148	560	565	91	208	21010	60738	-	-
1012	434	TRUE	TRUE	8141	4207	31577	73647	48	26	1154	1155	558	563	102	186	23436	69440	-	-
1013	475	TRUE	FALSE	6734	568	31909	72768	28	7	787	787	617	610	81	159	25175	72200	-	-
1014	1258	TRUE	FALSE	33867	3860	103057	193818	102	9	1052	1038	816	809	157	160	69190	189958	-	-
1015	1598	TRUE	FALSE	66130	4085	139638	243785	155	11	1052	1060	826	845	201	161	73508	239700	-	-
1016	194	TRUE	TRUE	5096	1073	19258	33083	51	11	1174	1170	733	730	124	176	14162	32010	-	-
1017	263	TRUE	TRUE	4164	2847	20207	45979	57	28	1177	1171	734	746	118	192	16043	43132	-	-
1018	581	TRUE	TRUE	9580	2190	31077	86435	58	11	1164	1165	479	479	95	156	21497	84245	-	-
1019	605	TRUE	TRUE	9991	3222	36006	95787	50	15	1225	1232	559	552	93	168	26015	92565	-	-
1020	460	TRUE	TRUE	10367	2131	27387	73431	65	8	1226	1223	563	563	102	163	17020	71300	-	-
1021	800	TRUE	TRUE	9593	3674	49593	126874	47	15	1189	1196	554	551	97	169	40000	123200	-	-
1022	675	TRUE	TRUE	12927	4357	46002	111682	61	17	1189	1194	560	562	110	176	33075	107325	-	-
1023	142	FALSE	TRUE	340	1156	8292	24160	5	15	1146	1144	594	594	61	177	7952	23004	-	-

Images:

A:1 Gel

B:2 Blot

Appendix - Audit Trail

time	user	entry
10/05/2018 10:08:29	Michael Hutchins	add image 1 Gel
10/05/2018 10:08:30	Michael Hutchins	add image 2 Blot
10/05/2018 10:09:09	Michael Hutchins	spots updated- new count: 795
10/05/2018 10:09:09	Michael Hutchins	detection sensitivity: 263
10/05/2018 10:09:09	Michael Hutchins	detection smoothing: 1
10/05/2018 10:09:09	Michael Hutchins	detection ROI: ROI Rectangle: 272-258 1110x1232
10/05/2018 10:09:09	Michael Hutchins	set base image: 1 Gel
10/05/2018 10:14:02	Michael Hutchins	perform alignment- 2 Blot to base image 1 Gel - 17 vectors
10/05/2018 10:14:38	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:14:38	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:14:42	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:14:42	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:14:44	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:14:44	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:14:52	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:14:52	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:15:09	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:15:09	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:15:26	Michael Hutchins	edit: delete 3 spot(s)
10/05/2018 10:15:42	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:15:46	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:15:50	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:15:54	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:15:58	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:16:00	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:16:00	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:16:02	Michael Hutchins	edit: delete 3 spot(s)
10/05/2018 10:16:02	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:16:08	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:16:08	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:16:10	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:16:10	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:16:19	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:16:22	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:16:25	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:16:28	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:16:32	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:16:32	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:16:34	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:16:34	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:16:38	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:16:38	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:16:41	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:16:41	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:16:47	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:16:47	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:16:53	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:16:53	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:16:57	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:16:57	Michael Hutchins	edit: add 1 spot(s)

10/05/2018 10:16:58	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:16:58	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:17:01	Michael Hutchins	edit: delete 4 spot(s)
10/05/2018 10:17:01	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:17:02	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:02	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:17:03	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:03	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:17:04	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:04	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:17:06	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:17:06	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:09	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:17:09	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:17:18	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:18	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:20	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:20	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:23	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:23	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:25	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:25	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:32	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:32	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:34	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:34	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:36	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:36	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:37	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:37	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:42	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:17:42	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:17:49	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:17:53	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:17:58	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:18:01	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:18:11	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:18:13	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:18:16	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:18:23	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:18:23	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:18:29	Michael Hutchins	edit: delete 4 spot(s)
10/05/2018 10:18:29	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:18:34	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:18:34	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:18:57	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:18:57	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:19:01	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:19:01	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:19:03	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:19:03	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:19:06	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:19:06	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:19:18	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:19:18	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:19:21	Michael Hutchins	edit: delete 3 spot(s)
10/05/2018 10:19:21	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:20:06	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:20:10	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:20:14	Michael Hutchins	edit: add 1 spot(s)

10/05/2018 10:20:16	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:20:21	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:20:28	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:20:41	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:20:41	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:20:43	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:20:43	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:20:48	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:20:48	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:20:54	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:20:57	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:21:08	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:21:08	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:21:11	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:21:11	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:21:14	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:21:14	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:21:19	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:21:39	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:21:59	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:22:01	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:22:04	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:22:07	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:22:09	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:22:17	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:22:17	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:22:18	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:22:18	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:22:23	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:22:23	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:22:25	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:22:25	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:22:32	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:22:32	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:22:35	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:22:35	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:22:41	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:22:41	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:22:53	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:22:55	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:23:03	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:23:05	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:23:13	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:23:13	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:23:14	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:23:14	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:23:17	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:23:17	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:23:20	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:23:20	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:23:39	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:23:41	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:23:45	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:23:45	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:23:52	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:24:15	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:24:16	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:24:23	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:24:23	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:24:27	Michael Hutchins	edit: add 1 spot(s)

10/05/2018 10:24:38	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:24:38	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:24:45	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:24:45	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:24:47	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:24:47	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:24:48	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:24:48	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:24:58	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:24:58	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:25:06	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:25:06	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:25:07	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:25:07	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:25:15	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:25:15	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:25:24	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:25:24	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:25:28	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:25:28	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:25:35	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:25:45	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:25:45	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:25:59	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:25:59	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:26:05	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:26:07	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:26:23	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:26:23	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:26:37	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:26:37	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:26:39	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:26:39	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:26:41	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:26:41	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:26:44	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:26:44	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:26:51	Michael Hutchins	edit: delete 3 spot(s)
10/05/2018 10:26:51	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:26:53	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:26:53	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:26:55	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:26:55	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:27:00	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:27:00	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:27:06	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:27:06	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:27:08	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:27:08	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:27:19	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:27:19	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:27:25	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:27:25	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:27:26	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:27:26	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:27:49	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:27:49	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:28:10	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:28:10	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:28:20	Michael Hutchins	edit: add 1 spot(s)

10/05/2018 10:28:27	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:28:27	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:28:35	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:28:43	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:28:44	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:29:20	Michael Hutchins	edit: delete 4 spot(s)
10/05/2018 10:29:24	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:29:26	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:29:29	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:29:30	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:29:33	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:29:40	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:29:40	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:29:43	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:29:47	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:29:53	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:29:55	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:29:59	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:30:04	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:30:07	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:30:07	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:30:16	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:30:20	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:30:21	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:30:23	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:30:23	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:30:29	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:30:38	Michael Hutchins	edit: delete 3 spot(s)
10/05/2018 10:30:41	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:31:58	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:31:58	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:32:05	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:32:08	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:32:15	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:32:15	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:32:39	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:33:09	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:33:20	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:33:27	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:33:33	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:34:49	Michael Hutchins	set 429 spots absent
10/05/2018 10:55:43	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:55:43	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:56:04	Michael Hutchins	edit: delete 3 spot(s)
10/05/2018 10:56:04	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:56:21	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:56:21	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:56:27	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 10:56:27	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 10:56:59	Michael Hutchins	set 9 spots absent
10/05/2018 10:57:04	Michael Hutchins	set 14 spots absent
10/05/2018 10:57:06	Michael Hutchins	set 17 spots absent
10/05/2018 10:57:07	Michael Hutchins	set 34 spots absent
10/05/2018 10:57:10	Michael Hutchins	set 4 spots absent
10/05/2018 10:57:12	Michael Hutchins	set 5 spots absent
10/05/2018 10:57:16	Michael Hutchins	set 19 spots absent
10/05/2018 10:57:17	Michael Hutchins	set 2 spots absent
10/05/2018 10:57:23	Michael Hutchins	set 9 spots absent
10/05/2018 10:57:25	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:30	Michael Hutchins	set 9 spots absent

10/05/2018 10:57:34	Michael Hutchins	set 20 spots absent
10/05/2018 10:57:36	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:36	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:37	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:39	Michael Hutchins	set 4 spots absent
10/05/2018 10:57:40	Michael Hutchins	set 5 spots absent
10/05/2018 10:57:41	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:41	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:42	Michael Hutchins	set 4 spots absent
10/05/2018 10:57:45	Michael Hutchins	set 7 spots absent
10/05/2018 10:57:47	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:47	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:50	Michael Hutchins	set 2 spots absent
10/05/2018 10:57:50	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:55	Michael Hutchins	set 23 spots absent
10/05/2018 10:57:56	Michael Hutchins	set 1 spots absent
10/05/2018 10:57:58	Michael Hutchins	set 8 spots absent
10/05/2018 10:57:59	Michael Hutchins	set 1 spots absent
10/05/2018 10:58:02	Michael Hutchins	set 18 spots absent
10/05/2018 10:58:05	Michael Hutchins	set 12 spots absent
10/05/2018 10:58:06	Michael Hutchins	set 2 spots absent
10/05/2018 10:58:08	Michael Hutchins	set 5 spots absent
10/05/2018 10:58:09	Michael Hutchins	set 5 spots absent
10/05/2018 10:58:14	Michael Hutchins	set 10 spots absent
10/05/2018 10:58:15	Michael Hutchins	set 7 spots absent
10/05/2018 10:58:17	Michael Hutchins	set 17 spots absent
10/05/2018 10:58:17	Michael Hutchins	set 2 spots absent
10/05/2018 10:58:20	Michael Hutchins	set 4 spots absent
10/05/2018 10:58:34	Michael Hutchins	set 19 spots absent
10/05/2018 10:58:36	Michael Hutchins	set 2 spots absent
10/05/2018 10:58:45	Michael Hutchins	set 5 spots absent
10/05/2018 10:58:48	Michael Hutchins	set 2 spots absent
10/05/2018 10:58:50	Michael Hutchins	set 1 spots absent
10/05/2018 10:58:55	Michael Hutchins	set 8 spots absent
10/05/2018 10:58:59	Michael Hutchins	set 4 spots absent
10/05/2018 10:59:00	Michael Hutchins	set 4 spots absent
10/05/2018 10:59:02	Michael Hutchins	set 4 spots absent
10/05/2018 10:59:03	Michael Hutchins	set 10 spots absent
10/05/2018 10:59:05	Michael Hutchins	set 4 spots absent
10/05/2018 10:59:06	Michael Hutchins	set 2 spots absent
10/05/2018 10:59:09	Michael Hutchins	set 7 spots absent
10/05/2018 10:59:11	Michael Hutchins	set 2 spots absent
10/05/2018 10:59:14	Michael Hutchins	set 2 spots absent
10/05/2018 10:59:19	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:21	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:21	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:22	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:23	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:23	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:26	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:27	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:27	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:32	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:37	Michael Hutchins	edit: delete 2 spot(s)
10/05/2018 10:59:37	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 10:59:41	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:42	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:42	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:43	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:43	Michael Hutchins	set 1 spots absent

10/05/2018 10:59:44	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:44	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:47	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:49	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:52	Michael Hutchins	set 1 spots absent
10/05/2018 10:59:56	Michael Hutchins	set 1 spots absent
10/05/2018 11:00:22	Michael Hutchins	set 4 spots absent
10/05/2018 11:00:23	Michael Hutchins	set 6 spots absent
10/05/2018 11:00:24	Michael Hutchins	set 8 spots absent
10/05/2018 11:00:27	Michael Hutchins	set 4 spots absent
10/05/2018 11:00:30	Michael Hutchins	set 12 spots absent
10/05/2018 11:01:00	Michael Hutchins	set 10 spots absent
10/05/2018 11:01:01	Michael Hutchins	set 1 spots absent
10/05/2018 11:01:04	Michael Hutchins	set 4 spots absent
10/05/2018 11:01:05	Michael Hutchins	set 2 spots absent
10/05/2018 11:01:05	Michael Hutchins	set 2 spots absent
10/05/2018 11:01:07	Michael Hutchins	set 6 spots absent
10/05/2018 11:01:10	Michael Hutchins	set 3 spots absent
10/05/2018 11:01:11	Michael Hutchins	set 1 spots absent
10/05/2018 11:01:13	Michael Hutchins	set 1 spots absent
10/05/2018 11:01:14	Michael Hutchins	set 1 spots absent
10/05/2018 11:01:16	Michael Hutchins	set 2 spots absent
10/05/2018 11:01:17	Michael Hutchins	set 1 spots absent
10/05/2018 11:01:18	Michael Hutchins	set 9 spots absent
10/05/2018 11:01:21	Michael Hutchins	set 18 spots absent
10/05/2018 11:01:25	Michael Hutchins	set 1 spots present
10/05/2018 11:01:28	Michael Hutchins	set 1 spots absent
10/05/2018 11:01:29	Michael Hutchins	set 2 spots absent
10/05/2018 11:01:39	Michael Hutchins	set 2 spots absent
10/05/2018 11:01:43	Michael Hutchins	set 5 spots absent
10/05/2018 11:01:44	Michael Hutchins	set 4 spots absent
10/05/2018 11:01:49	Michael Hutchins	set 1 spots absent
10/05/2018 11:01:50	Michael Hutchins	set 1 spots absent
10/05/2018 11:01:51	Michael Hutchins	set 1 spots absent
10/05/2018 11:01:53	Michael Hutchins	set 5 spots absent
10/05/2018 11:02:18	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:16	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:17	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:23	Michael Hutchins	set 8 spots absent
10/05/2018 11:03:25	Michael Hutchins	set 4 spots absent
10/05/2018 11:03:28	Michael Hutchins	set 2 spots absent
10/05/2018 11:03:38	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:38	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:38	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:48	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:49	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:52	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:53	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:53	Michael Hutchins	set 1 spots absent
10/05/2018 11:03:54	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:03	Michael Hutchins	set 3 spots absent
10/05/2018 11:04:24	Michael Hutchins	set 14 spots absent
10/05/2018 11:04:26	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:29	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:31	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:31	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:32	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:34	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:34	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:34	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:36	Michael Hutchins	set 1 spots absent
10/05/2018 11:04:37	Michael Hutchins	set 1 spots absent

10/05/2018 11:04:43	Michael Hutchins	set 13 spots absent
10/05/2018 11:04:52	Michael Hutchins	set 4 spots absent
10/05/2018 11:04:57	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 11:04:57	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 11:05:09	Michael Hutchins	edit: delete 1 spot(s)
10/05/2018 11:05:09	Michael Hutchins	edit: add 2 spot(s)
10/05/2018 11:05:14	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:15	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:17	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:20	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:22	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:26	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:40	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:41	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:43	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:45	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:47	Michael Hutchins	set 1 spots absent
10/05/2018 11:05:48	Michael Hutchins	set 1 spots absent
10/05/2018 11:06:20	Michael Hutchins	edit: add 1 spot(s)
10/05/2018 11:06:30	Michael Hutchins	set 2 spots present
10/05/2018 11:06:32	Michael Hutchins	set 2 spots present
10/05/2018 11:06:38	Michael Hutchins	set 6 spots present
10/05/2018 11:06:40	Michael Hutchins	set 2 spots present
10/05/2018 11:07:11	Michael Hutchins	set 3 spots present